

Service Manual for:

NL Millennium **Series**

Public Use Wheelchair Lifts

Series C5

DOT — Public Use Lift

"DOT — Public Use Lift" verifies that this platform lift meets the "public use lift" requirements of FMVSS No. 403. This lift may be installed on all vehicles appropriate for the size and weight of the lift, but must be installed on buses, school buses, and multi-purpose passenger vehicles other than motor homes with a gross vehicle weight rating (GVWR) that exceeds 4,536 kg (10,000 lb).




"Providing Access to the World"®

International Corporate Hdqrs: P.O. Box 310 Winamac, IN 46996 USA
1-800-THE LIFT® (574) 946-6153 FAX: (574) 946-4670



! WARNING

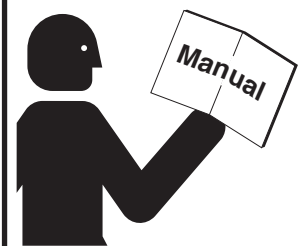

Read manual
before installing
or servicing lift.
Failure to do so
may result in
serious bodily
injury and/or
property damage.

Braun NL Series

33995 Rev. C

December

2009

Congratulations

We at The Braun Corporation wish to express our fullest appreciation on your new purchase. With you in mind, our skilled craftsmen have designed and assembled the finest lift available.

This manual provides service-related material. Refer to the FMVSS No. 403 Quick Reference Installation Sheet for installation instructions, operating instructions and maintenance procedures.

Braun Millennium Series™ lifts are built for dependability and will provide years of pleasure and independence as long as the lift is installed and serviced as specified by a Braun certified technician, and the lift is operated by an instructed person.

Sincerely,
THE BRAUN CORPORATION




Ralph W. Braun
Chief Executive Officer

Warranty and Registration Instructions

Immediately upon receiving the lift, examine the unit for any damage. Notify the carrier at once with any claims.


Two warranty/registration cards (shown right) are located in the lift-mounted manual storage pouch. The sales representative must process one of the cards. The consumer must fill out the other card and mail it to The Braun Corporation. The warranty is provided on the back cover of this manual. The **warranty cards must be processed to activate the warranty.**

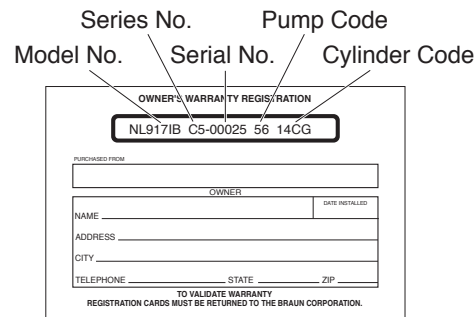

The card is titled "OWNER'S WARRANTY REGISTRATION". It contains a box with the following information: NL9171B C5-00025 56 14CG. Below this, there are fields for "PURCHASED FROM", "OWNER", "NAME", "ADDRESS", "CITY", "TELEPHONE", "STATE", and "ZIP". There is also a "DATE INSTALLED" field. At the bottom, it says "TO VALIDATE WARRANTY REGISTRATION CARDS MUST BE RETURNED TO THE BRAUN CORPORATION."

Sample Warranty/Registration Card

Two Braun Serial No./Series No. identification tags (shown below) are posted on the lift. One I.D. tag is posted on the opposite pump side vertical arm. A second I.D. tag is located on the opposite pump side tower. Both I.D. tags provide the product identification information provided on the warranty/registration card. Record the information in the space provided (or document on a copy). **This information must be provided when filing a warranty claim or ordering parts.**

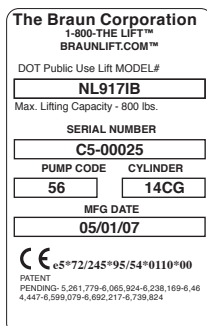

The tag is from "The Braun Corporation" with phone number "1-800-THE LIFT™" and website "BRAUNLIFT.COM™". It identifies the "DOT Public Use Lift MODEL# NL9171B" with a "Max. Lifting Capacity - 800 lbs.". It lists the "SERIAL NUMBER C5-00025", "PUMP CODE 56", and "CYLINDER 14CG". The "MFG DATE" is "05/01/07". It also includes a CE mark and patent information: "eS*72/245*95/54*0110*00", "PATENT PENDING- 5,261,779-6,065,924-6,238,169-6,46 4,447-6,599,079-6,692,217-6,739,824".

Sample Serial No./Series No. Identification Tag

Model No. _____

Series No. _____

Serial No. _____

Pump Code _____

Cylinder Code _____

Date of Manufacture _____

Contents

Troubleshooting and Maintenance

Lift Terminology	2
Switch and Sensor Locations	3
Certification Checklist Diagnostic Procedures	4
Platform Fold Pressure Adjustment	5
Outer Barrier Fold Pressure Adjustment	5
Platform Angle Adjustment	6-7
Platform Stop Blocks	7
Platform Floor Level Adjustment	8
Bridging Microswitch Adjustment	9
LCD Lift Codes	10-11
Lubrication Diagram	12
Maintenance and Lubrication Schedule	13-16
Lift Electrical Schematic	17
Lift Wiring Diagram	18

Hydraulics

Hydraulic Schematic	19
Hydraulics Parts List	20
Hydraulics Diagram	21

Repair Parts

Pump Module

Pump Module Parts List	22
Pump Module Diagram	23

Lift Exploded Views

42" Floor-to-Ground Lift Models

(NL915IB, NL915FIB, NL918IB, and NL918FIB)

Repair Parts List	24
Exploded View.....	25

48" Floor-to-Ground Lift Models


(NL916IB, NL916FIB, NL917IB, NL917FIB,
NL919IB and NL919FIB)

Exploded View	26
Repair Parts List	27


Warranty

Braun® Limited Warranty	28-30
-------------------------------	-------

Lift Terminology

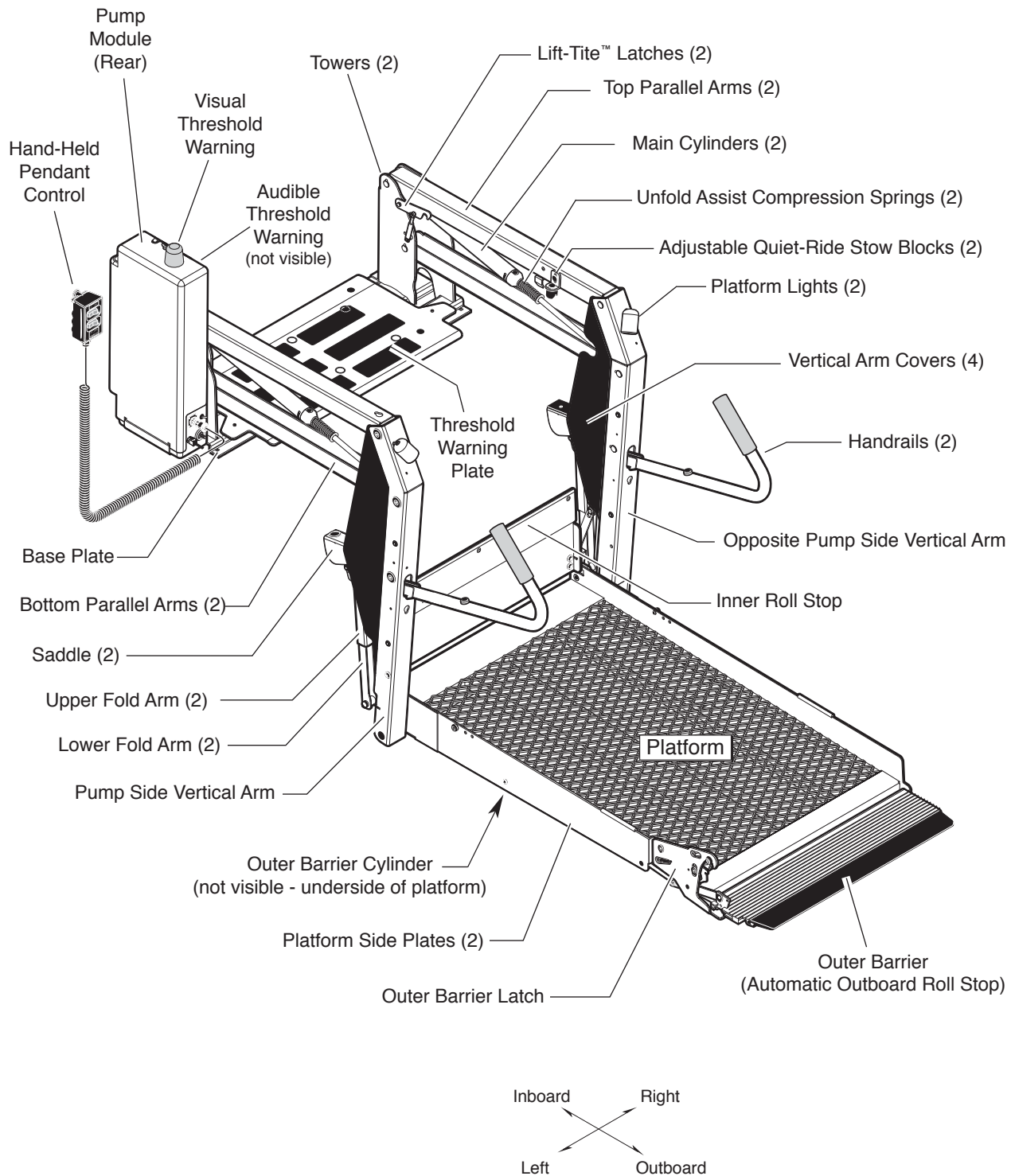

Switch and Sensor Locations

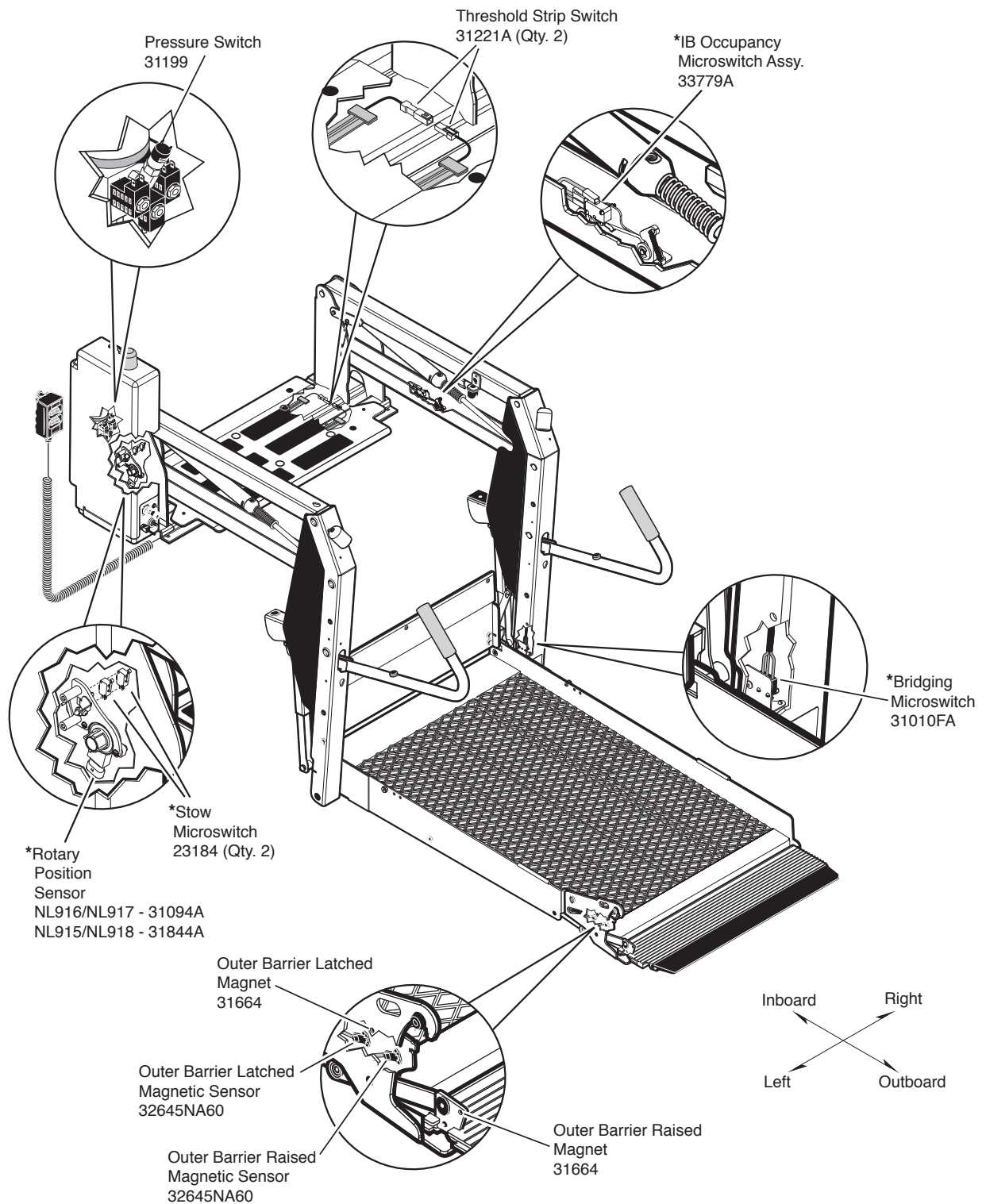

*Note: Mirror image for right (front) pump lifts.


Certification Checklist Diagnostic Procedures

The following operations and conditions must be functionally verified in order for the lift to be FMVSS 403/404 compliant. If an operation does not function as described or a condition is not met, follow the referenced procedures to correct the problem or contact a Braun Corporation Product Support representative at 1-800-THE LIFT®.

- Vehicle movement is prevented unless the lift door is closed, ensuring the lift is stowed.
 1. Verify on the pump module mounted interlock connector that the lift not stowed signal (pin 5) has an open signal - OR - lift stowed signal (pin 7) has a +12 volt signal - OR - (pin 9) has a ground signal (depends on interlock used).
 2. Refer to the interlock installation instructions.
- Lift operation shall be prevented unless the vehicle is stopped and vehicle movement is prevented.
 1. Verify vehicle secure signal (pin 6) has a +12 volt signal.
 2. Refer to the interlock installation instructions.
- The platform will not fold/stow if occupied.
 - Refer to Platform Fold Pressure Adjustment procedures.
- The inner roll stop will not raise if occupied.
 - Call Product Support
- The outer barrier will not raise if occupied.
 - Refer to Outer Barrier Fold Pressure Adjustment procedures.
- Verify platform lighting when lift is deployed and pendant illumination when lift is powered.
 1. Replace bulb(s) in the light housing.
 2. Check fuse (5 amp fuse on circuit board; F13)
- A visual and audible warning will activate if the threshold area is occupied when the platform is at least 1" below floor level.
 1. Remove the threshold warning plate
 2. Verify the threshold strip switch connectors are connected
 3. Replace the threshold strip switch
 4. Reinstall the threshold warning plate
- Platform movement is prohibited beyond the position where the inner roll stop is fully deployed (up).
 - Call Product Support.
- Lift platform movement shall be interrupted unless the outer barrier is deployed (up).
 - Call Product Support.


Platform Fold Pressure Adjustment

1. Position the platform at the floor level loading position.
2. Loosen the hex nut on the adjustment screw (do not remove hex nut).
3. Turn the adjustment screw counter clockwise until the platform does not fold when the Fold button is pressed.
4. Turn the adjustment screw clockwise in 1/4 turn increments and press the Fold button until the platform folds completely (Note: Return the platform to floor level position after each attempt to fold the platform).
5. Turn the adjustment screw an additional 1/8 turn after the platform folds successfully.
6. Tighten the hex nut without moving the adjustment screw.
7. Verify the platform will not stow while occupied.

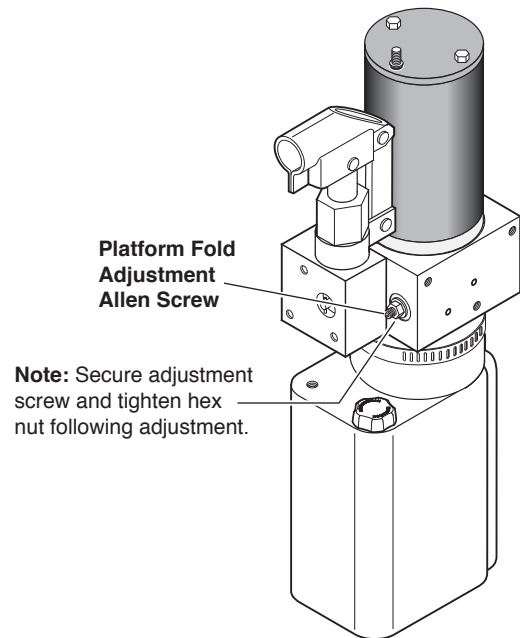

Outer Barrier Fold Pressure Adjustment

1. Lower the platform to the ground level loading position.
2. Remove the pump cover.
3. Turn the adjustment screw (shown below) counter clockwise until the outer barrier does not raise when the Up button is pressed.
4. Turn the adjustment screw clockwise in 1/4 turn increments until the outer barrier raises and fully locks in position when the Up button is pressed. DO NOT turn the adjustment screw beyond this position.
5. Verify the outer barrier will not raise when occupied.
6. Reattach the pump cover.

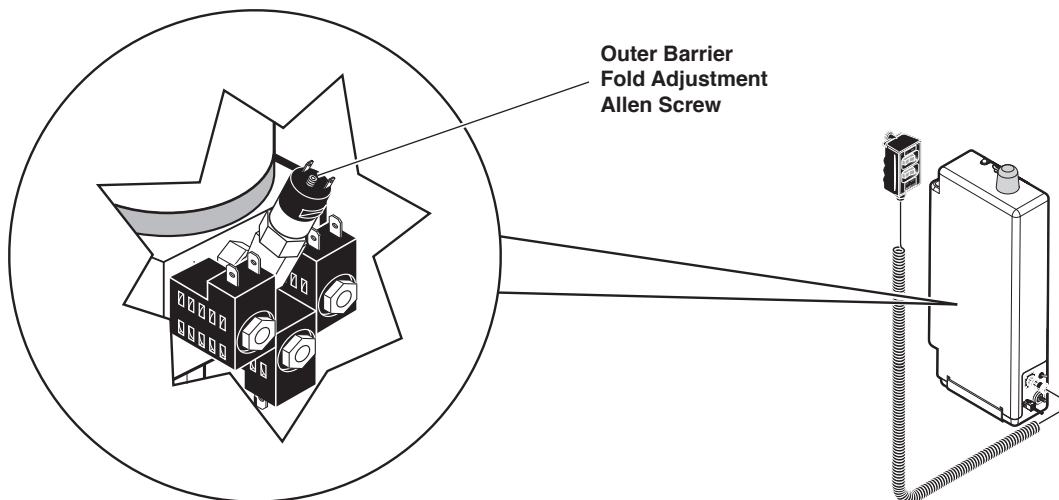

Platform Angle Adjustment

Adjusting the platform angle based on the relationship of the platform at ground level directly affects the angle of the platform when positioned at floor level.

Unfold the lift and visually examine the angle of the platform when positioned at floor level. Lower the platform fully and note the angle of the platform when it reaches ground level also. The **outboard end of the platform (toe) must contact the ground first** when the platform is lowered (all lift models).

Millennium “NL” Series:

The platform angle **should** be adjusted so there is a balance between the angle at both positions (**equal amount of angle**). **Angle A** should **equal angle B** as shown in Figure C.

Century “NCL” Series:

Vista “NVL” Series:

The platform angle **must** be adjusted so the outboard end of the platform (toe) is angled down slightly when positioned at floor level. See Figure D. The outboard end of the platform **must** contact the ground first to **ensure** the spring-loaded outer **barrier unfolds fully**.

Adjustment Procedures: Adjustment Allen screws are provided on each side of the lift platform for adjusting the platform angle. Details and photos are provided on the opposite page.


Base Plate Wedges: Installations where base plate wedges are used require more platform angle adjustment than normal.

Platform Stop Blocks: When adjusting platform angle, ensure both **stop blocks are making full contact with the vertical arms** (details on opposite page).

Floor Level Adjustment: Following platform angle adjustment, set floor level positioning as detailed in **Platform Floor Level Adjustment** (details on page 8).

Millennium “NL” Series


Figure C

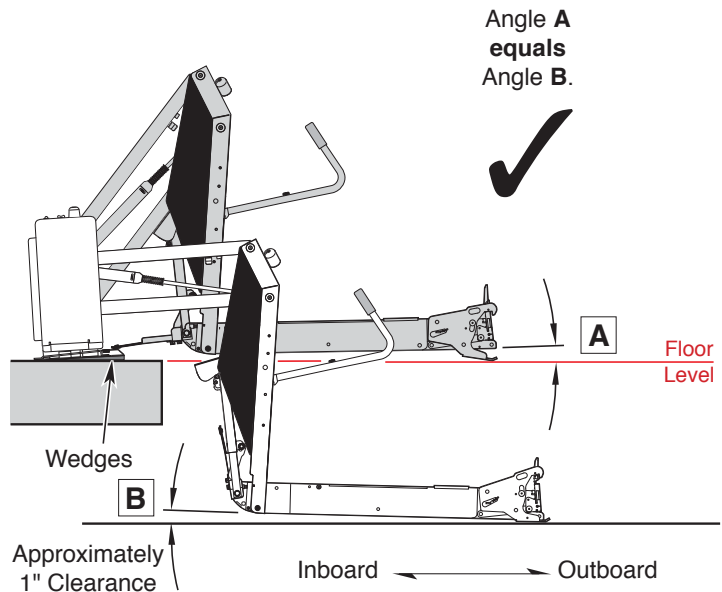

Century “NCL” Series

Vista “NVL” Series

Figure D

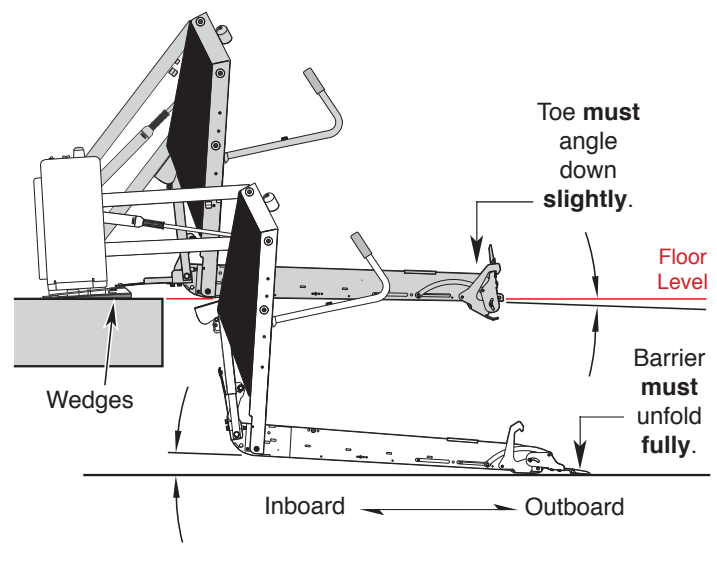

Platform Angle Adjustment


Adjustment Allen screws are provided on each side of the lift platform for adjusting the platform angle. Adjust platform angle as specified on previous page.

To **raise** the **outboard end** of platform - turn adjustment screw **clockwise**.


To **lower** the **outboard end** of platform - turn adjustment screw **counterclockwise**.

Note: Both adjustment screws **must be adjusted equally**.

Apply Loctite® to adjustment screws following adjustment.

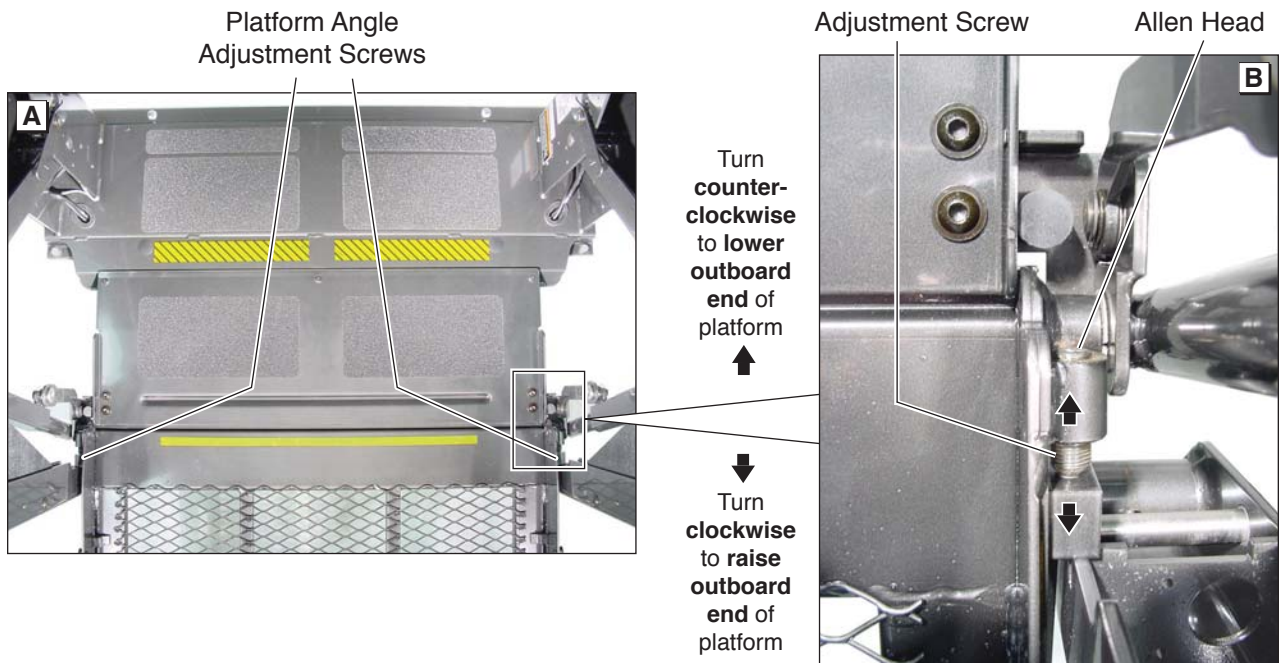

Platform Stop Blocks

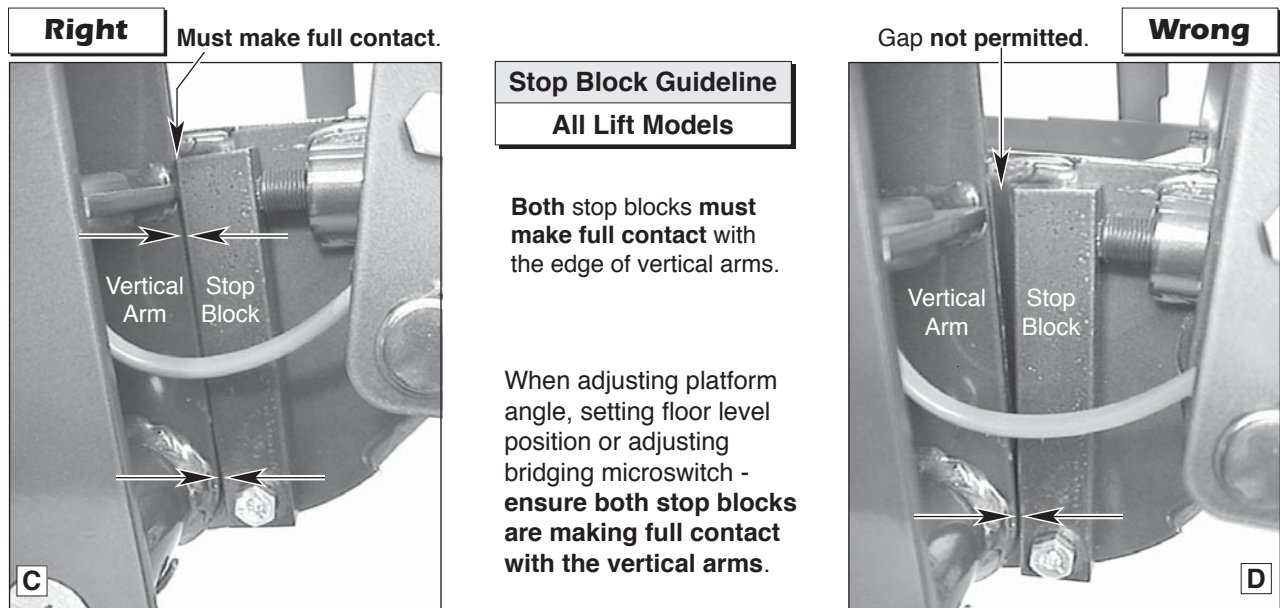

Platform Floor Level Adjustment

Before setting floor level position:

- Adjust platform angle as detailed on page 6.
- Ensure both **stop blocks are making full contact with vertical arms** (details on page 7).

Floor Level Adjustment:

1. Position platform at desired floor level position (passenger loading/unloading height).

Position platform such that:

- a. the inner roll stop is laying flat on the threshold plate
- b. platform has **not** begun to fold

Note: Use hand pump to position platform at proper posi-

tion if unable to stop platform when powering lift.

2. Turn Lift Power switch Off.
3. Press Floor Position Set button (located between pump housing and lift tower).
4. While pressing the Floor Position Set button, turn the Lift Power switch On.
5. Continue pressing the Floor Position Set button until the lift sounds three “beeps.”
6. Release the Floor Position Set button.
7. Cycle lift to verify that platform stops at the set floor level position. **Note:** If platform does not stop at



Floor Level Set Button

the intended position - repeat adjustment procedures. If repeating procedures fails - refer to Diagnostics section below.

Diagnostics: Diagnostic codes have been established in event the lift platform floor position does not set (the lift does not sound three “beeps” - see Step 5 above). The control board located inside the pump hous-

ing is equipped with an LCD screen. Remove the pump cover to access the LCD screen. The following diagnostic codes will help resolve floor position setting problems:

Century/Vista Jumper



- 91 – The platform position is out of a predetermined acceptable range
- 92 – The Bridging Microswitch is not activated (adjust switch or lower the platform)
- 93 – The Inner Roll Stop Occupied switch is not activated (adjust switch)
- 94 – The Outer Barrier Up switch is not activated (adjust switch)
- 95 – The Outer Barrier Latched sensor is not activated (Century and Vista: Verify jumper is installed on the outboard barrier latch switch – see photo at right; Millennium: Check latch)

LCD Screen





Bridging Microswitch Adjustment

Bridging Microswitch:

The bridging microswitch is located at the bottom of the right (front) vertical arm. See Photos A and B. An adjustment screw is built into the platform stop block.

The bridging microswitch will be deactivated if the outboard end of the platform contacts the unloading surface before the inboard end while the DOWN switch is pressed.

The deactivated bridging microswitch interrupts power to the DOWN switch circuit and stops the down travel of the platform (while allowing the outer barrier to deploy).

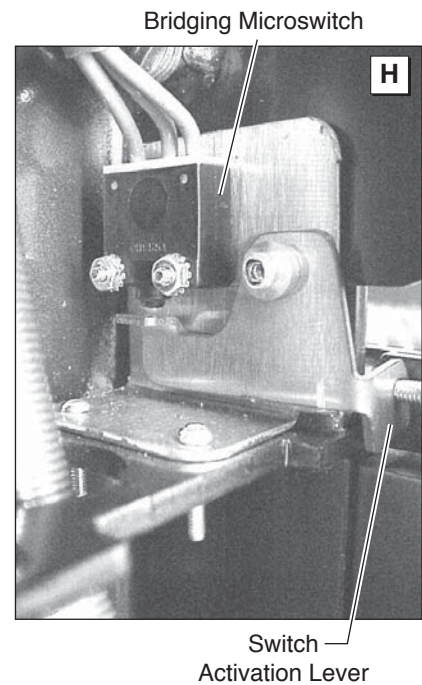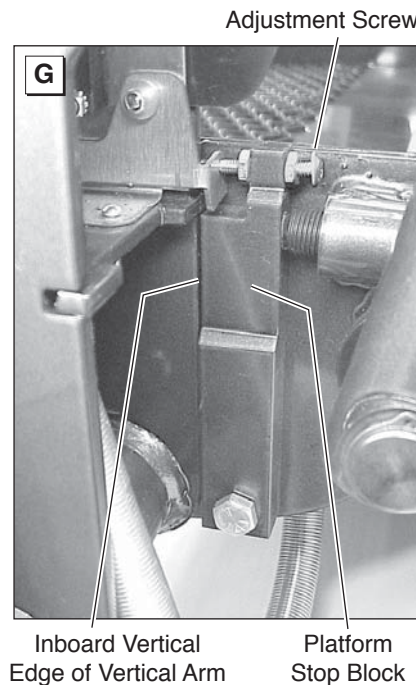

Adjustment Guidelines:

Caution: This adjustment is factory set and typically should not require adjustment.

Before adjusting the bridging microswitch:

- Adjust platform angle as detailed on page 6.

- Adjust floor level position as detailed on page 8.
- **Ensure** platform is in the fully unfolded position.
- **Ensure** both stop blocks are making full contact with vertical arms (details on page 7).

Failure to follow these guidelines will result in an overtight adjustment that will bend the adjustment screw, break the bridging switch and/or result in lift operation failure.

Bridging Microswitch Adjustment

1. Position the platform at floor level (or slightly below). Ensure the platform is not partially folded. From the stowed position, press the UNFOLD switch until the pump stops running and the platform comes to a full stop. From below floor level, press the UP switch until the pump stops running and the platform comes to a full stop.
2. Check **both** platform stop blocks to **ensure there is no gap** between the stop blocks and the vertical arms. See Photos B, C and D on page 7. A gap may indicate that the platform is not fully unfolded or the platform angle adjustment bolts are not adjusted equally on both sides.
3. Loosen the two jam nuts on the bridging switch activation screw (see Photo G). Tighten the screw until the activation lever is fully depressed. Verify that there is no play in the activation lever and tighten the jam nuts.
4. Test lift for proper operation.

LCD Lift Codes

To better understand the Braun LCD Trouble shooting display you must first understand the numbers that appear on the screen. There are Flashing Codes, Solid Error Codes, and Solid Normal Operational Codes.

Flashing Codes #65-89: About 10 seconds after an operation has stopped there are a set of scrolling flashing numbers that indicate whenever a particular sensor or switch has been activated. These numbers will start at number 65 and scroll up to number 89, then start the sequence over. Remember they are not error codes. Keep in mind that the lift will display codes for different positions and certain flashing codes must be present for that position, you will not always get an error code.

Solid Error Codes #50-64: These are the numbers that will come on the screen when the audible and visual alarm goes off, and will direct you to where the problem exists. These numbers will only stay on the screen for about 10 seconds and then the flashing codes will scroll indicating what sensors are active. This sequence of codes will keep repeating. It is important to be looking at the screen when trying to get the lift to fail.

Solid Normal Operational Codes: There are also solid numbers that will appear while and after the lift is moving that indicate the lift operation and platform position.

Troubleshooting Procedures

1. While **looking** at the LCD screen, operate the lift until the failure occurs. Read the number that comes on the display the moment the alarm goes off and the light starts to flash. This code will only stay on the screen for 10 seconds.
2. Look up the number on the correct error code sheet and determine what part on the lift is causing the failure. Go to the part on the lift that is suspected of causing the failure and look for anything obvious like magnets missing, broken wires etc. If nothing is found, the next step is to determine if that sensor is sending a signal to the board.
3. Bring the platform to the level that the sensor should be activated using the backup pump if needed. At this point, look up the flashing code that corresponds to that sensor in the error code sheet, look at the LCD screen and wait for the flashing scrolling numbers to appear. If the number is not included in the scrolling numbers, you know that sensor is the problem. You should then check the harness or try another magnet with the south side of the magnet facing the sensor and see if the number will come up on the display.
4. If the problem is still not found or the harness is suspected, the voltages should to be checked to and from the sensor to find the exact location of the problem. First determine the 3 wire colors for this sensor at the board and understand the 3 voltage readings needed to operate the sensor, the 12V power, 8V power, and the 8V input signal to the PC board when activated by a magnet. First check for a 8V input signal coming from the sensor to the wire going into the PC board, if there is 11V on this wire, the sensor is not being activated by the magnet. Next check the 12V and 8V wires at the PC board plug leading to the sensor. Once verified at the plug on the PC board, the voltages should then be checked at the next plug down on the harness going to the sensor until the location of the problem is found.

Anytime you see the code for that specific switch you will have 8 volts on that colored wire on the 8 or 18 pin connector from that switch. IE: Outboard Barrier is closed "72" will appear on the screen and also 8V will be present on the signal wire from that switch, if no code is present the voltage will be 11V.

LCD Lift Codes

Listed below are codes that the lift controller outputs during lift operation. The codes will be displayed on an LCD screen located on the lift control board inside the pump module. See the Manual Operating Instructions in the operator's manual for pump cover removal instructions.

Non-Flashing Numbers


- 01 – Platform stowed
- 02 – Platform unfolding
- 03 – Platform unfolding paused
- 04 – Platform at floor level
- 05 – Platform beginning to lower
- 06 – Platform lowering (threshold cannot be occupied from this point down)
- 07 – Outer barrier moving to horizontal position
- 08 – Platform at ground level
- 09 – Outer barrier moving to vertical position
- 10 – Platform raising
- 11 – Platform raising paused at floor
- 12 – Platform folding (limited pressure)
- 13 – Platform folding (full pressure)
- 14 – Timed fold (cinching lift tite) or (anti-rattle state)
- 15 – Platform folding stopped
- 16 – Paused fold
- 17 – Platform between ground and 3" above ground
- 18 – Platform above 3"
- 19 – Outer barrier moving to horizontal position
- 23 – Outer barrier going back down after occupant detected
- 28 – Illegal function/not defined
- 29 – Interlock fault not recognized (or has been cleared but a motion button is still pressed)
- 30 – Platform location unknown
- 31 – Platform location transition state; attempting to locate position
- 35 – Two or more motion buttons are being pressed
- 36 – The retention belt cannot be buckled while trying to fold or unfold
- 37 – Motion button being pressed is not a valid motion
- 50 – Outer barrier is not up above inboard barrier locked position
- 51 – Threshold is occupied when platform is 1" or more below floor level
- 52 – Inner roll stop is not up and locked below inner roll stop locked position
- 53 – Inner roll stop occupied sensor is not activated between floor and inner roll stop up position
- 54 – Outer barrier is occupied before it is up
- 55 – Outer barrier is not latched when above the inner roll stop locked position (Millennium only)
- 56 – Outer barrier is not up and latched and bridge switch did not deactivate
- 57 – Outer barrier is not up and latched and ground detect switch did not deactivate (Century and Vista only)
- 58 – Outer barrier is not up and latched and the platform is 3" above the ground
- 59 – Outer barrier is not up after pausing platform travel
- 60 – Verify the IB occupied switch is functioning correctly when lift is below floor level. The kickout gas springs might be worn, replace before using.
- 75 – Low voltage detected; must turn off power switch to reset LCD
- 77 – Vehicle secure interlock has not been activated
- 90 – Position will be set if you keep holding the button until it beeps
- 91 – Position is out of a predetermined acceptable range of floor position
- 92 – Bridge switch is not made, needs adjusting
- 93 – Inner roll stop occupied switch is not made, position needs to be moved or switch should be adjusted
- 94 – Outer barrier is not made, fix and try again
- 95 – Outer barrier latch is not made (check for jumper on Century and Vista lifts, check latch on Millennium lifts)
- 99 – Controller program is not valid; replace controller

Flashing Numbers

- Flashing 65 – Unfold button is pressed
- Flashing 66 – Fold button is pressed
- Flashing 67 – Down button is pressed
- Flashing 68 – Up button is pressed
- Flashing 69 – Bridge switch is activated
- Flashing 70 – Outer barrier latch switch is activated
- Flashing 71 – Ground detect switch is activated
- Flashing 72 – Outer barrier up switch is activated
- Flashing 73 – Inner roll stop up switch is activated
- Flashing 74 – Inner roll stop occupied switch is activated
- Flashing 76 – Outboard barrier occupied switch is activated
- Flashing 78 – Threshold tape switch "A" is activated
- Flashing 79 – Threshold tape switch "B" is activated
- Flashing 80 – Position set button is pressed

Maintenance and Lubrication

Lubrication Diagram

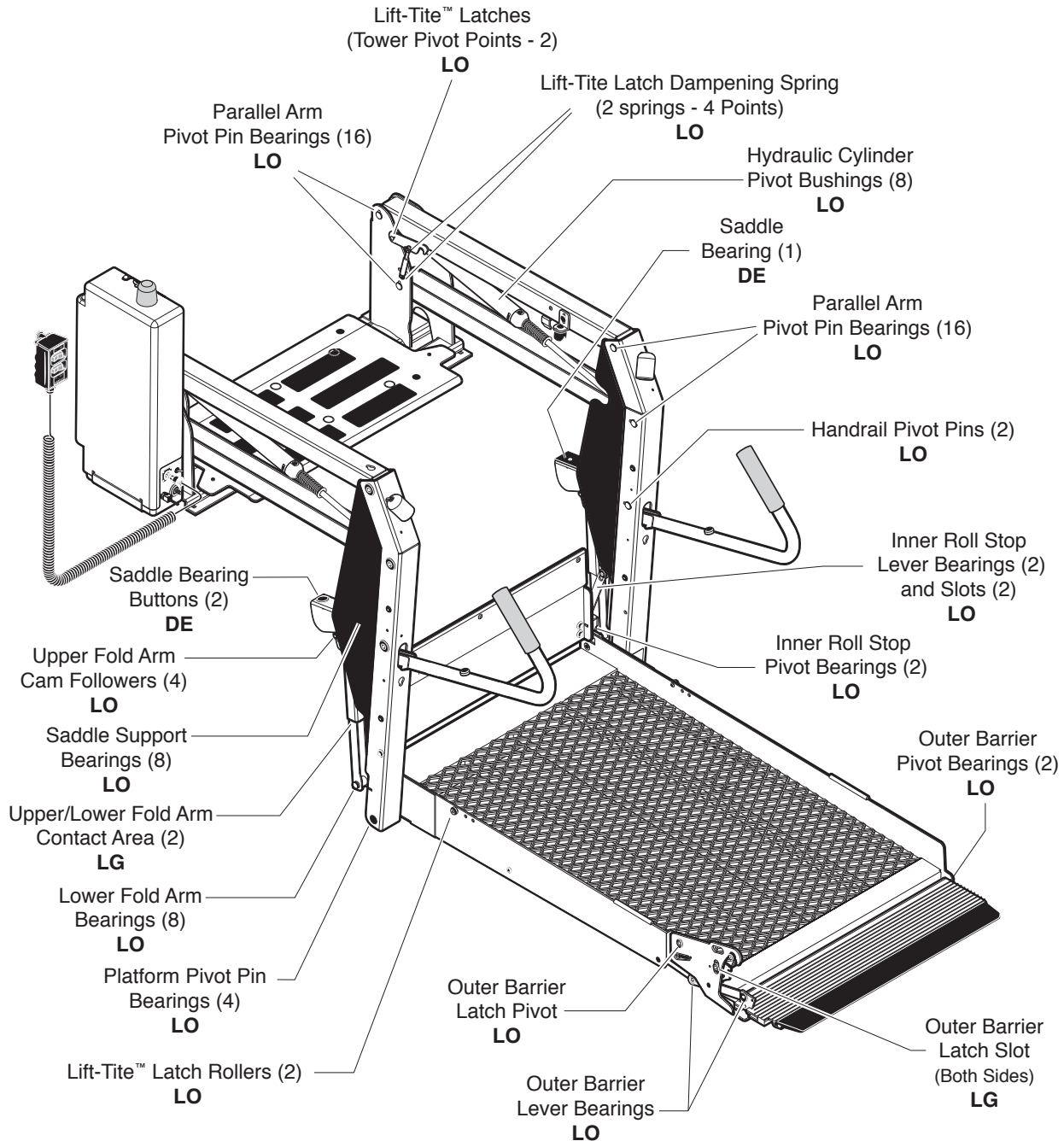

See the Maintenance/Lubrication Schedule for recommended applications per number of cycles.

Lubricant	Type	Specified (recommended) Lubricant	Available Amount	Braun Part No.
LO - Light Oil	Light Penetrating Oil (30 weight or equivalent)	LPS2, General Purpose Penetrating Oil	11 oz. Aerosol Can	15807
DE - Door-Ease	Stainless Stick Style (tube)	Door-Ease Stick (tube)	1.68 oz.	15806
LG - Light Grease	Light Grease (Multipurpose)	Lubriplate	14 oz. Can	15805

Maintenance and Lubrication Schedule

Proper maintenance is necessary to ensure safe, troublefree operation. Inspecting the lift for any wear, damage or other abnormal conditions should be a part of all transit agencies's daily service program. Simple inspections can detect potential problems.

The maintenance and lubrication procedures specified in this schedule **must** be performed by a Braun authorized service representative at the scheduled intervals according to the number of cycles.

Braun dual parallel arm lifts are equipped with hardened pins and self-lubricating bushings to decrease wear, provide smooth operation and extend the service life of the lift.

When servicing the lift at the recommended intervals, inspection and lubrication procedures specified in the previous sections should be repeated. **Clean** the components and the surrounding area **before applying lubricants**. LPS2 General Purpose Penetrating Oil is recommended where Light Oil is called out. Use of improper lubricants can attract dirt or other contaminants which could result in wear or damage to the components. Platform components exposed to contaminants when lowered to the ground may require extra attention.

Lift components requiring grease are lubricated during assembly procedures. When these components are replaced, grease must be applied during installation procedures. Specified lubricants are available from The Braun Corporation (part numbers provided above).

All listed inspection, lubrication and maintenance procedures should be repeated at "750 cycle" intervals fol-

lowing the scheduled "4500 Cycles" maintenance. These intervals are a general guideline for scheduling maintenance procedures and will vary according to lift use and conditions. Lifts exposed to severe conditions (weather, environment, contamination, heavy usage, etc.) may require inspection and maintenance procedures to be performed more often than specified.

Maintenance Indicator: The Lift Ready green LED mounted on top of the pump cover will begin to blink after every 750 cycles. The blinking LED will not affect the functions of the lift, but is a reminder to complete necessary maintenance and lubrication.

Once the lift has been serviced, fully stow the lift. Once stowed, press the UP button on the hand pendant and the Floor Level Set button on the back side of the pump cover until the Lift Ready green LED stops blinking.

Discontinue lift use immediately if maintenance and lubrication procedures are not properly performed, or if there is any sign of wear, damage or improper operation. Contact your sales representative or call The Braun Corporation at 1-800-THE LIFT®. One of our national Product Support representatives will direct you to an authorized service technician who will inspect your lift.

WARNING

Maintenance and lubrication procedures must be performed as specified by an authorized service technician. Failure to do so may result in serious bodily injury and/or property damage.

750 Cycles	Outer barrier pivot bearings (2)	Apply Light Oil - See Lubrication Diagram
	Outer barrier latch pivot point	Apply Light Oil - See Lubrication Diagram
	Outer barrier latch slot	Apply Light Grease to both sides of slot. See Lubrication diagram
	Outer barrier lever bearings (2)	Apply Light Oil - See Lubrication Diagram
	Lift-Tite™ latches (tower pivot points - 2)	Apply Light Oil - See Lubrication Diagram
	Lift-Tite™ latch gas (dampening) spring pivot points (2 springs - 4 points)	Apply Light Oil - See Lubrication Diagram
	Inspect Lift-Tite™ latches and gas springs for wear or damage (bent, deformed or misaligned), positive securement (external snap rings) and proper operation	Resecure, replace damaged parts or otherwise correct as needed. Note: Apply Light Grease to Lift-Tite™ latch tower pivot point if replacing latch.
	Inspect outer barrier for proper operation	Correct or replace damaged parts.
continued	Inspect outer barrier latch for proper operation, positive securement, and detached or missing spring	Correct or replace damaged parts and/or relubricate. See Lubrication Diagram

Maintenance and Lubrication Schedule

<div>continued</div> <div>750 Cycles</div>	<table> <tr> <td>Adjust Fold Pressure And Outer Barrier Fold Pressure</td><td>See Platform Fold Pressure Adjustment and Outer Barrier Fold Pressure Adjustment</td></tr> <tr> <td>Verify FMVSS 403/404 Certification Checklist</td><td>See Certification Checklist Diagnostic Procedures</td></tr> <tr> <td>Inspect lift for wear, damage or any abnormal condition</td><td>Correct as needed.</td></tr> <tr> <td>Inspect lift for rattles</td><td>Correct as needed.</td></tr> </table>	Adjust Fold Pressure And Outer Barrier Fold Pressure	See Platform Fold Pressure Adjustment and Outer Barrier Fold Pressure Adjustment	Verify FMVSS 403/404 Certification Checklist	See Certification Checklist Diagnostic Procedures	Inspect lift for wear, damage or any abnormal condition	Correct as needed.	Inspect lift for rattles	Correct as needed.																								
Adjust Fold Pressure And Outer Barrier Fold Pressure	See Platform Fold Pressure Adjustment and Outer Barrier Fold Pressure Adjustment																																
Verify FMVSS 403/404 Certification Checklist	See Certification Checklist Diagnostic Procedures																																
Inspect lift for wear, damage or any abnormal condition	Correct as needed.																																
Inspect lift for rattles	Correct as needed.																																
<div>1500 Cycles</div> <div>continued</div>	<div>Perform all procedures listed in previous section also</div> <table> <tr> <td>Upper/lower fold arms (2)</td><td>Apply grease (synthetic) to contact areas between upper/lower fold arms. See Lubrication Diagram</td></tr> <tr> <td>Platform pivot pin bearings (4)</td><td>Apply Light Oil - See Lubrication Diagram</td></tr> <tr> <td>Lower fold arm bearings (8)</td><td>Apply Light Oil - See Lubrication Diagram</td></tr> <tr> <td>Inner roll stop pivot bearings (2)</td><td>Apply Light Oil - See Lubrication Diagram</td></tr> <tr> <td>Inner roll stop lever bearings (2)</td><td>Apply Light Oil - See Lubrication Diagram</td></tr> <tr> <td>Inner roll stop lever slot (2)</td><td>Apply Light Oil - See Lubrication Diagram</td></tr> <tr> <td>Saddle support bearings (8)</td><td>Apply Light Oil - See Lubrication Diagram</td></tr> <tr> <td>Upper fold arm cam followers (4)</td><td>Apply Light Oil - See Lubrication Diagram</td></tr> <tr> <td>Parallel arm pivot pin bearings (16)</td><td>Apply Light Oil - See Lubrication Diagram</td></tr> <tr> <td>Handrail pivot pins (2)</td><td>Apply Light Oil - See Lubrication Diagram</td></tr> <tr> <td>Hydraulic cylinder pivot bushings (8)</td><td>Apply Light Oil - See Lubrication Diagram</td></tr> <tr> <td>Inspect Lift-Tite™ latch rollers for wear or damage, positive securement and proper operation (2)</td><td>Correct, replace damaged parts and/or relubricate.</td></tr> <tr> <td>Inspect inner roll stop for: • Wear or damage • Proper operation. Roll stop should just rest on top surface of the threshold plate. • Positive securement (both ends)</td><td>Resecure, replace or correct as needed. See Platform Angle Instructions and Platform Floor Level Adjustment.</td></tr> <tr> <td>Inspect handrail components for wear or damage, and for proper operation</td><td>Replace damaged parts.</td></tr> <tr> <td>Inspect microswitches for securement and proper adjustment.</td><td>Resecure, replace or adjust as needed. See Microswitch Adjustment Instructions.</td></tr> <tr> <td>Make sure lift operates smoothly</td><td>Realign towers and vertical arms. Lubricate or correct as needed.</td></tr> </table>	Upper/lower fold arms (2)	Apply grease (synthetic) to contact areas between upper/lower fold arms. See Lubrication Diagram	Platform pivot pin bearings (4)	Apply Light Oil - See Lubrication Diagram	Lower fold arm bearings (8)	Apply Light Oil - See Lubrication Diagram	Inner roll stop pivot bearings (2)	Apply Light Oil - See Lubrication Diagram	Inner roll stop lever bearings (2)	Apply Light Oil - See Lubrication Diagram	Inner roll stop lever slot (2)	Apply Light Oil - See Lubrication Diagram	Saddle support bearings (8)	Apply Light Oil - See Lubrication Diagram	Upper fold arm cam followers (4)	Apply Light Oil - See Lubrication Diagram	Parallel arm pivot pin bearings (16)	Apply Light Oil - See Lubrication Diagram	Handrail pivot pins (2)	Apply Light Oil - See Lubrication Diagram	Hydraulic cylinder pivot bushings (8)	Apply Light Oil - See Lubrication Diagram	Inspect Lift-Tite™ latch rollers for wear or damage, positive securement and proper operation (2)	Correct, replace damaged parts and/or relubricate.	Inspect inner roll stop for: • Wear or damage • Proper operation. Roll stop should just rest on top surface of the threshold plate. • Positive securement (both ends)	Resecure, replace or correct as needed. See Platform Angle Instructions and Platform Floor Level Adjustment.	Inspect handrail components for wear or damage, and for proper operation	Replace damaged parts.	Inspect microswitches for securement and proper adjustment.	Resecure, replace or adjust as needed. See Microswitch Adjustment Instructions.	Make sure lift operates smoothly	Realign towers and vertical arms. Lubricate or correct as needed.
Upper/lower fold arms (2)	Apply grease (synthetic) to contact areas between upper/lower fold arms. See Lubrication Diagram																																
Platform pivot pin bearings (4)	Apply Light Oil - See Lubrication Diagram																																
Lower fold arm bearings (8)	Apply Light Oil - See Lubrication Diagram																																
Inner roll stop pivot bearings (2)	Apply Light Oil - See Lubrication Diagram																																
Inner roll stop lever bearings (2)	Apply Light Oil - See Lubrication Diagram																																
Inner roll stop lever slot (2)	Apply Light Oil - See Lubrication Diagram																																
Saddle support bearings (8)	Apply Light Oil - See Lubrication Diagram																																
Upper fold arm cam followers (4)	Apply Light Oil - See Lubrication Diagram																																
Parallel arm pivot pin bearings (16)	Apply Light Oil - See Lubrication Diagram																																
Handrail pivot pins (2)	Apply Light Oil - See Lubrication Diagram																																
Hydraulic cylinder pivot bushings (8)	Apply Light Oil - See Lubrication Diagram																																
Inspect Lift-Tite™ latch rollers for wear or damage, positive securement and proper operation (2)	Correct, replace damaged parts and/or relubricate.																																
Inspect inner roll stop for: • Wear or damage • Proper operation. Roll stop should just rest on top surface of the threshold plate. • Positive securement (both ends)	Resecure, replace or correct as needed. See Platform Angle Instructions and Platform Floor Level Adjustment.																																
Inspect handrail components for wear or damage, and for proper operation	Replace damaged parts.																																
Inspect microswitches for securement and proper adjustment.	Resecure, replace or adjust as needed. See Microswitch Adjustment Instructions.																																
Make sure lift operates smoothly	Realign towers and vertical arms. Lubricate or correct as needed.																																


Maintenance and Lubrication Schedule

<p>continued</p> <p>1500 Cycles</p>	<table> <tr> <td>Inspect external snap rings: • Lower fold arm (6) • Lift-Tite™ latch (2) • Lift-Tite™ latch roller (2) • Lift-Tite™ latch gas (dampening) spring (4) • Upper fold arm cam followers (4) • Outer barrier hydraulic cylinder mounting pin (2) • Inner roll stop lever bracket pins (2)</td><td>Resecure or replace if needed.</td></tr> <tr> <td>Inspect inner roll stop locks (2) and torson springs (2) for wear or damage and for proper operation.</td><td>Replace damaged parts. Apply Light Oil to inner roll stop lock pivot point.</td></tr> <tr> <td>Inspect lower fold arm pins (2), axles (2) and bearings (8) for wear or damage and positive securement</td><td>Replace damaged parts and resecure as needed. Apply Light Oil.</td></tr> <tr> <td>Remove pump module cover and inspect: • Hydraulic hoses, fittings and connections for wear or leaks • Harness cables, wires, terminals and connections for securement or damage • Control board, circuit breaker, power switch and lights for securement or damage</td><td>Resecure, replace or correct as needed.</td></tr> </table>	Inspect external snap rings: • Lower fold arm (6) • Lift-Tite™ latch (2) • Lift-Tite™ latch roller (2) • Lift-Tite™ latch gas (dampening) spring (4) • Upper fold arm cam followers (4) • Outer barrier hydraulic cylinder mounting pin (2) • Inner roll stop lever bracket pins (2)	Resecure or replace if needed.	Inspect inner roll stop locks (2) and torson springs (2) for wear or damage and for proper operation.	Replace damaged parts. Apply Light Oil to inner roll stop lock pivot point.	Inspect lower fold arm pins (2), axles (2) and bearings (8) for wear or damage and positive securement	Replace damaged parts and resecure as needed. Apply Light Oil.	Remove pump module cover and inspect: • Hydraulic hoses, fittings and connections for wear or leaks • Harness cables, wires, terminals and connections for securement or damage • Control board, circuit breaker, power switch and lights for securement or damage	Resecure, replace or correct as needed.												
Inspect external snap rings: • Lower fold arm (6) • Lift-Tite™ latch (2) • Lift-Tite™ latch roller (2) • Lift-Tite™ latch gas (dampening) spring (4) • Upper fold arm cam followers (4) • Outer barrier hydraulic cylinder mounting pin (2) • Inner roll stop lever bracket pins (2)	Resecure or replace if needed.																				
Inspect inner roll stop locks (2) and torson springs (2) for wear or damage and for proper operation.	Replace damaged parts. Apply Light Oil to inner roll stop lock pivot point.																				
Inspect lower fold arm pins (2), axles (2) and bearings (8) for wear or damage and positive securement	Replace damaged parts and resecure as needed. Apply Light Oil.																				
Remove pump module cover and inspect: • Hydraulic hoses, fittings and connections for wear or leaks • Harness cables, wires, terminals and connections for securement or damage • Control board, circuit breaker, power switch and lights for securement or damage	Resecure, replace or correct as needed.																				
<p>4500 Cycles</p> <p>continued</p>	<table> <tr> <td colspan="2">Perform all procedures listed in previous section also</td></tr> <tr> <td>Inspect cotter pins on platform pivot pin (2)</td><td>Resecure, replace or correct as needed</td></tr> <tr> <td>Hydraulic Fluid (Pump) - Check level. Note: Fluid should be changed if there is visible contamination. Inspect the hydraulic system (cylinder, hoses, fittings, seals, etc.) for leaks if fluid level is low.</td><td>Use Braun 32840-QT (Exxon® Univis HVI 26) hydraulic fluid (do not mix with Dextron III or other hydraulic fluids). Check fluid level with platform lowered fully and roll stop unfolded fully. Fill to within 1/2" of the bottom of the 1 1/2" fill tube (neck).</td></tr> <tr> <td>Inspect cylinders, fittings and hydraulic connections for wear, damage or leaks</td><td>Tighten, repair or replace if needed.</td></tr> <tr> <td>Inspect outer barrier cylinder hose assembly (hose, fasteners, connections, etc.) for wear, damage or leakage</td><td>Tighten, repair or replace if needed.</td></tr> <tr> <td>Inspect parallel arms, bearings and pivot pins for visible wear or damage</td><td>Replace if needed.</td></tr> <tr> <td>Inspect parallel arm pivot pin mounting bolts (8)</td><td>Tighten or replace if needed.</td></tr> <tr> <td>Inspect platform pivot pins, bearings and vertical arms for wear, damage and positive securement</td><td>Replace damaged parts and resecure as needed. Apply Light Grease during reassembly procedures.</td></tr> <tr> <td>Inspect upper/lower fold arms, saddle, saddle support and associated pivot pins and bearings for visible wear or damage</td><td>Replace if needed.</td></tr> <tr> <td>Inspect gas springs (cylinders - 4) for wear or damage, proper operation and positive securement</td><td>Tighten, replace or correct as needed</td></tr> </table>	Perform all procedures listed in previous section also		Inspect cotter pins on platform pivot pin (2)	Resecure, replace or correct as needed	Hydraulic Fluid (Pump) - Check level. Note: Fluid should be changed if there is visible contamination. Inspect the hydraulic system (cylinder, hoses, fittings, seals, etc.) for leaks if fluid level is low.	Use Braun 32840-QT (Exxon® Univis HVI 26) hydraulic fluid (do not mix with Dextron III or other hydraulic fluids). Check fluid level with platform lowered fully and roll stop unfolded fully . Fill to within 1/2" of the bottom of the 1 1/2" fill tube (neck).	Inspect cylinders, fittings and hydraulic connections for wear, damage or leaks	Tighten, repair or replace if needed.	Inspect outer barrier cylinder hose assembly (hose, fasteners, connections, etc.) for wear, damage or leakage	Tighten, repair or replace if needed.	Inspect parallel arms, bearings and pivot pins for visible wear or damage	Replace if needed.	Inspect parallel arm pivot pin mounting bolts (8)	Tighten or replace if needed.	Inspect platform pivot pins, bearings and vertical arms for wear, damage and positive securement	Replace damaged parts and resecure as needed. Apply Light Grease during reassembly procedures.	Inspect upper/lower fold arms, saddle, saddle support and associated pivot pins and bearings for visible wear or damage	Replace if needed.	Inspect gas springs (cylinders - 4) for wear or damage, proper operation and positive securement	Tighten, replace or correct as needed
Perform all procedures listed in previous section also																					
Inspect cotter pins on platform pivot pin (2)	Resecure, replace or correct as needed																				
Hydraulic Fluid (Pump) - Check level. Note: Fluid should be changed if there is visible contamination. Inspect the hydraulic system (cylinder, hoses, fittings, seals, etc.) for leaks if fluid level is low.	Use Braun 32840-QT (Exxon® Univis HVI 26) hydraulic fluid (do not mix with Dextron III or other hydraulic fluids). Check fluid level with platform lowered fully and roll stop unfolded fully . Fill to within 1/2" of the bottom of the 1 1/2" fill tube (neck).																				
Inspect cylinders, fittings and hydraulic connections for wear, damage or leaks	Tighten, repair or replace if needed.																				
Inspect outer barrier cylinder hose assembly (hose, fasteners, connections, etc.) for wear, damage or leakage	Tighten, repair or replace if needed.																				
Inspect parallel arms, bearings and pivot pins for visible wear or damage	Replace if needed.																				
Inspect parallel arm pivot pin mounting bolts (8)	Tighten or replace if needed.																				
Inspect platform pivot pins, bearings and vertical arms for wear, damage and positive securement	Replace damaged parts and resecure as needed. Apply Light Grease during reassembly procedures.																				
Inspect upper/lower fold arms, saddle, saddle support and associated pivot pins and bearings for visible wear or damage	Replace if needed.																				
Inspect gas springs (cylinders - 4) for wear or damage, proper operation and positive securement	Tighten, replace or correct as needed																				

Maintenance and Lubrication Schedule

<div>continued</div> <div>4500 Cycles</div>	<div>Perform all procedures listed in previous section also</div> <div> <div>Inspect saddle bearing (1) and saddle bearing buttons (2)</div> <div>Apply Door-Ease or replace if needed. See Lubrication Diagram.</div> </div> <div> <div>Inspect vertical arm plastic covers</div> <div>Resecure or replace if needed.</div> </div> <div> <div>Inspect power cable</div> <div>Resecure, repair or replace if needed.</div> </div> <div> <div>Mounting</div> <div>Check to see that the lift is securely anchored to the vehicle and there are no loose bolts, broken welds, or stress fractures.</div> </div> <div> <div>Decals and Antiskid</div> <div>Replace decals if worn, missing or illegible. Replace antiskid if worn or missing.</div> </div>
<div>Consecutive 750 Cycle Intervals</div>	<div>Repeat all previously listed inspection, lubrication and maintenance procedures at 750 cycle intervals.</div>

Lift Electrical Schematic

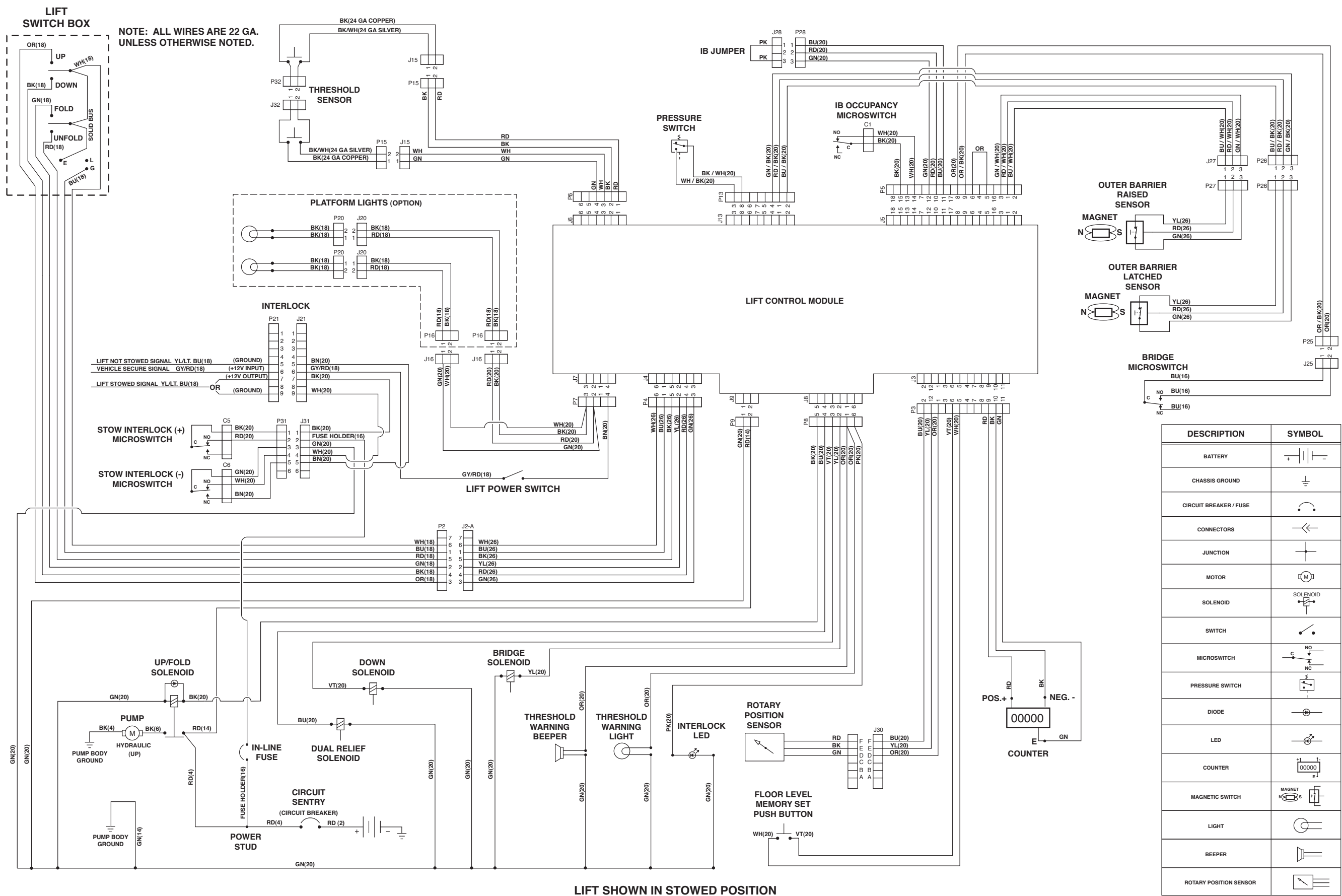


LIFT SHOWN IN STOWED POSITION

**Unfold for:
Lift Electrical
Schematic**



Hydraulic Schematic

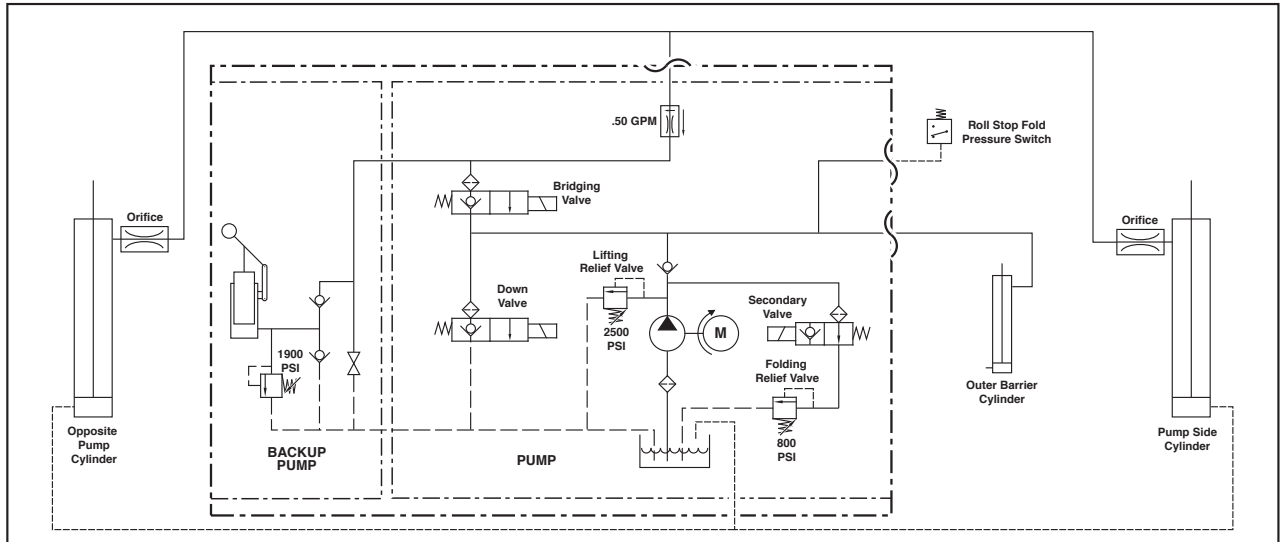

Description	Symbol	Description	Symbol
Fixed Displacement Pump		Hydraulic Port	
Pump Motor		2 Way 2 Position Solenoid Valve	
Backup Pump		Pressure Compensated Flow Control	
Single Acting Cylinder		Relief Valve	
Check Valve		Filter Screen	
Unfold Orifice		Vented Reservoir	
Manual Shutoff Valve		Orifice	
Pressure Switch			

Hydraulics Parts List


Part Numbers of Items Dedicated per Lift Model (Front and Rear Pump Models)						
Item	Qty.	Description	NL915	NL916	NL917	NL918
1	1	Hose Assembly, 1/8" (Opposite Pump Side)	16004A-081	16004A-079	16004A-086	16004A-084
2	1	Hose Assembly, 1/8" (Pump Side)	16004A-046	16004A-046	16004A-046	16004A-046
3	2	Cylinder ✓ (C1514.3-9408 shown)	C1514.3-9407	C1514.3-9408	C1514.3-9408	C1514.3-9407
4	1	(Rear Pump) Hose Assembly, 3/16" w/Two Guards	915-2601-123	915-2601-139	915-2601-139	915-2601-123
5	1	(Front Pump) Hose Assembly, 3/16" w/Two Guards	915-5602-82	915-5602-92.5	915-5602-92.5	915-5602-82
6	1	(Front Pump) Hose Assembly, 3/16"	915-5603-80.5	915-5603-83.5	915-5603-86.5	915-5603-83.5
Part Numbers of Items Identical on all Lift Models (Front and Rear Pump Models)						
Item	Qty.	Description	Part No.			
7	1	Pump Assembly (M268-0113)12V - with Reservoir/with Backup Pump)	31038-12V			
8	1	Solenoid, Up - 4-Post - Trombetta (Flat Mount)	31374			
9	1	Power Cable, Up Solenoid to Motor	29049			
10	1	Motor, Pump - 12 Volt - Low RPM	31350			
11	1	Elbow, 65° - Male 7/16-20 SAE O-Ring to 7/16-20 SAE O-Ring Boss	26578			
12	2	Valve Assembly, "Down / Bridge" (complete)	31348K			
13	3	Coil (only), "Down / Bridge" Valve - (shown below)	31122			
14	2	Cartridge (only), "Down / Bridge" Valve - (shown below)	26078			
15	1	Valve Assembly, "Duel Relief" (complete)	31120K			
16	1	Cartridge (only), "Duel Relief" Valve - (shown below)	31121			
17	1	Cap, Reservoir Filler	16935			
18	4	O-Ring (only), Hand Pump Mounting	17351			
19	1	Hand Pump (Backup) with O-Rings (Item 18)	26074			
20	3	Screw, 1/4-20 x 2-1/4", Allen Head	26080			
21	1	Clamp, Reservoir - H-48 (M259)	17069			
22	1	Reservoir, Hydraulic Fluid	16934			
23	1	Fitting, 90° - 1/8" NPT x 1/8" Barb - Plastic	87563			
24	1	Switch, Pressure, 7/16-20 SAE O-Ring Male	31199			
25	2	Fitting, Male 7/16-20 SAE O-Ring to Male 7/16-20 JIC 37°	24504			
26	1	Elbow, (1) Female 7/16-20 JIC 37 Swivel - (2) Male 7/16-20 JIC 37°	26579			
27	2	Elbow, 90°, 7/16-20 SAE O-Ring Male - 7/16-20 JIC 37° Male, Orifice	26667			
28	2	Elbow, 90°, 1/4" NPT Male to 1/4" Barbed	15150			
29	1	Hose, Thermal Plastic - Black, 1/8" I.D.	23742R*			
30	1	Connector, Plastic "Y", 1/8" O.D.	18877			
31	1	Fitting, 90° - Male 9/16-18 SAE O-Ring to Male 7/16-20 JIC 37°	87622			
32	1	Cylinder, Roll Stop (915-3318NA Contains Items 31,32 & 35)	25309NA			
33	1	Handle with Grip	17206A			
34	1	Diode Assembly, Up Solenoid	73906A			
35	1	Breather Vent, 1/4" NPT	26200-2			
36	1	Kit, Hydraulic Port Service Cap	27049K			

✓ **Seal Kits:** If repairing a cylinder, order Seal Kit #1500-0500P.


* Raw material items ordered and priced per inch (order specified length).

Hydraulic Fluid
When adding or changing hydraulic fluid, use Braun 32840-QT (Exxon® Univis HVI 26) hydraulic fluid (do not mix with Dextron III or other hydraulic fluids).


12 "Down and Bridge" Valve (complete)

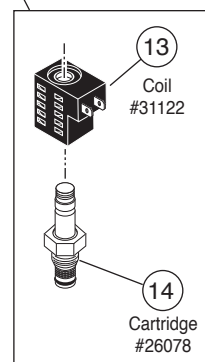

15 "Duel Relief" Valve (complete)

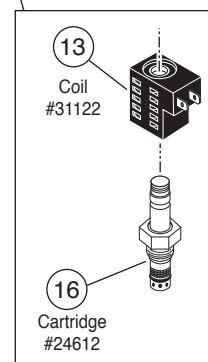

Hydraulics Diagram

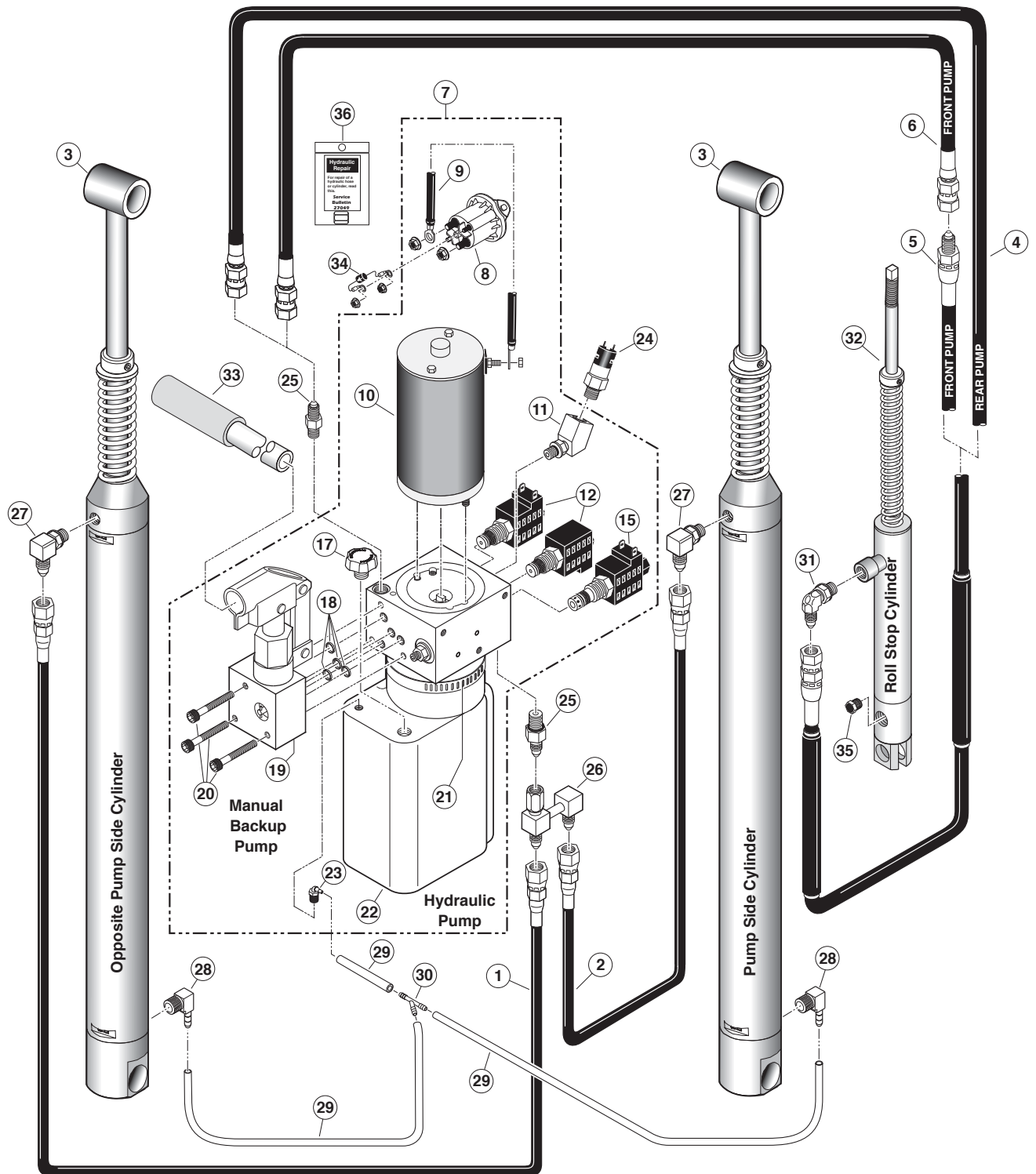

Pump Module Parts List

Part Numbers of Items Dedicated per Lift Model			NL915IB, NL916IB, NL917IB, NL918IB and NL919IB	NL915FIB, NL916FIB, NL917FIB, NL918FIB and NL919FIB
Item	Qty.	Description		
	1	Pump Module (complete), 12 Volt, Rear	917-5516RNA	917-5516FNA
1	1	Plate, Backing / Mounting - Rear	985-2501RN	985-2501FN
2	1	Cover, Pump, 2-Piece - Back / Bottom (915-0513RNAPT or 915-0513FNAPT Includes Items 2 and 16 - 22)	915-0513RNPT	915-0513FNPT
3	1	Cover, Pump - 2-Piece - Top / Front (915-0519RNAPT or 915-0519FNAPT Includes Items 3 and 39 - 52)	915-0519RNPT	915-0519FNPT
Part Numbers of Items Identical on all Lift Models				
Item	Qty.	Description		Part No.
4	1	Pump Assembly (M268-0113) 12V-120G - Dual Relief / ARS (Includes Items 5 & 6)		31038-12V
5	1	Power Cable, Up Solenoid to Motor		29049
6	1	Solenoid, Up - 4-Post - Trombetta (Flat Mount)		31374
7	1	Diode Assembly, Up Solenoid		73906A
8	2	Nut, #10-32, Serrated Flange		83080
9	2	Fitting, Male 7/16-20 O-Ring to Male 7/16-20 JIC 37°		24504
10	1	Elbow, Female Swivel 7/16-20 JIC 37° to (2) Male 7/16-20 JIC 37°		26579
11	1	Switch, Pressure		31199
12	1	Control, Hand Pendant Assembly - Electronic - Non-Shielded Cable - Coiled ♦See note below		33658A
13	1	Control Board Assembly		100159-001
14	8	Standoff, Snap-In		31011
15	1	Switch, Push Button		31753
16	1	Recepticle, Clip On		28803
17	1	Lens, Threshold Warning - Red		30704
18	1	Decal, Warning / Pressure Relief Valve (Not shown)		22249
19	1	Spacer, Lens - NHTSA		31386
20	1	Metal Ring Base - Lamp		30971
21	1	Socket, Lamp		30703
22	3	Screw, #8-32 x 1/2" Pan Head Phillips - Thread Cut		30974
23	1	Bulb, Light		19841
24	1	Diode, Green LED, Panel Mount		29545
25	4	Rivet, Snap - .124" Hole x .197"/.236" Thick		29874
26	2	Washer, 5/16" Flat		10063
27	2	Bolt, 5/16-18 x 3/4", Nylock, Hex * See note below		29608
28	2	Bolt, 5/16-18 x 1/2", Nylock, Hex * See note below		10012
29	1	Cable, Ground		22166A
30	1	Washer, 5/16" External Tooth		16368
31	1	Beeper, Constant - High Output		33251
32	1	Switch, Toggle		31787
33	1	Decal, Lift Power - On/Off (Not shown)		21494
34	1	Stud, Power Feed		26084
35	1	Rubber Boot, Red ♦See note below		82046
36	1	Harness, Lift Interlock Connection		32639A
37	1	Wire Assembly, Lift Stowed Connection ♦See note below		32638A
38	1	Pump Handle with Grip		17206A
39	3	Rivet, Pop, SD43BS - 1/8" - .13"/.19"		12954
40	1	Plug, Window - Clear		30443
41	1	Stud, Wing Head - 1/4 Turn		28804
42	1	Retainer, Push On		28805
43	1	Washer, Nylon, 1/4" I.D. x 11/16" O.D. x .030"		12690
44	1	Clip, Pump Handle - Top		915-5517
45	1	Clip, Pump Handle - Bottom		915-5518
46	1	Decal, Manual Instructions - Public (Not shown)		31412
47	1	Decal, Removal / Installaton - Pump Cover (Not shown)		29051
48	1	Decal, Warning - Control Board Damage - ESD (Not shown)		30787
49	1	Decal, Removal / Installation Pump Handle (Not shown)		29052
50	1	Decal, Installation & Warning - Hydraulic Pressure Switch (Not shown)		27154
51	1	Decal, Dual Relief Adjustment (Not shown)		32201
52	1	Decal, LCD Lift Codes (Not shown)		32410
53	1	Cable, Pump Module Power Hookup (Not shown - see Wiring Diagram)		26082A-4
54	1	Harness, Power - Stow Interlock (Not shown - see Wiring Diagram)		915-4530RNA
55	1	Harness, Interlock / Lighting (Not shown - see Wiring Diagram)		915-4533NA
56	1	Harness, Up / Down Solenoid (Not shown - see Wiring Diagram)		915-2533NA
57	1	Harness, Threshold Switch Extension (Not shown - see Wiring Diagram)		915-4534NA
58	1	Harness, Hand Control - Pump Housing to Circuit Board (Not shown - see Wiring Diagram)		100245-001


* Apply red #271 Thread Locker Loctite® to the four hex bolts (items 27 and 28) if a blue nylon patch is not present on the bolts when retrofitting an M268 pump assembly. Loctite® is available from The Braun Corporation under part number 11522-1.

◆ Indicates items available for replacement part purposes **only**. These items **are not** included with replacement pump modules.

Pump Module Diagram

Pump Mounting Bolts

Apply red #271 Thread Locker Locktite® to the three pump mounting bolts (items 27 and 28) if a blue nylon patch is not present on the bolts when retrofitting an M268 pump assembly. Locktite® is available from The Braun Corporation under part number 11522-1.

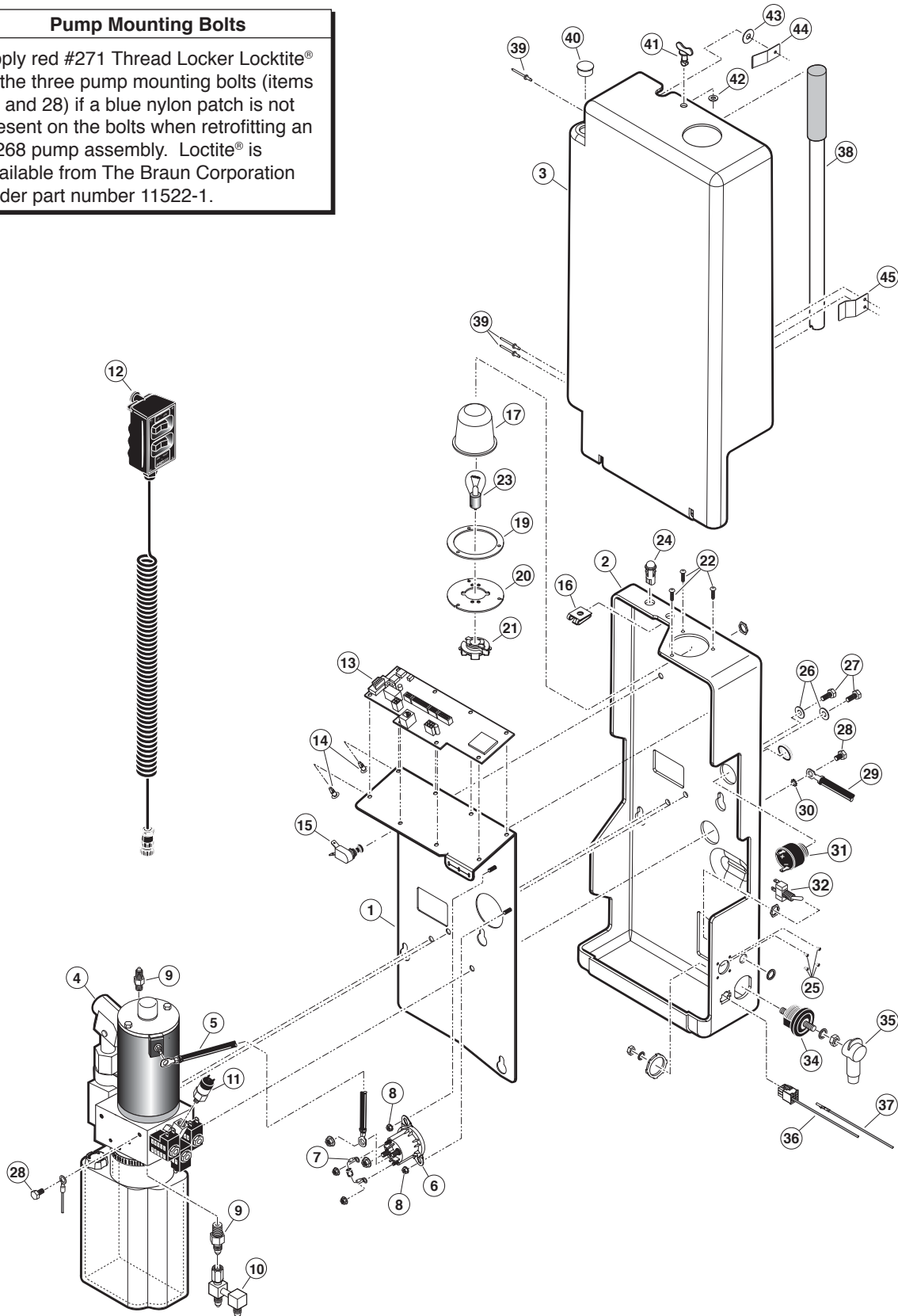


Note: Rear pump shown, front pump mirror image.

42" Floor-to-Ground Lift Repair Parts List


Part Numbers of Items Dedicated per Lift Model						
Item	Qty.	Description	NL915IB	NL915FIB	NL918IB	NL918FIB
1	1	Base Weldment	915R5142NW30	915F5142NW30	918R5142NW33	918F5142NW33
2	1	Cover, Pump, 2-Piece, Back-Bottom	915-0513RNPT	915-0513FNPT	915-0513RNPT	915-0513FNPT
3	1	Cover, Pump, 2-Piece, Top-Front Assy.	915-0519RNPT	915-0519FNPT	915-0519RNPT	915-0519FNPT
4	1	Platform Weldment	917-53050NW	917-53050NW	917-53350NW	917-53350NW
5	1	Outer Barrier Assembly (Inc. 5,14,22,51,168)	915-0312NA30Y	915-0312NA30Y	915-0312NA33Y	915-0312NA33Y
6	1	Inner Roll Stop	918-0147NWWY-30	918-0147NWWY-30	918-0147NWWY	918-0147NWWY
7	1	Bridge Switch Assembly	31010FA	31010FA	31010FA	31010FA
8	1	Plate, Backing/Mounting	985-2501RN	985-2501FN	985-2501RN	985-2501FN
9	1	Base Cover Assy., (Traction Tape Not Shown)	915-5142CNA30Y	915-5142CNA30Y	915-5142CNA33Y	915-5142CNA33Y
10	1	Parallel Arm Assy., Bottom, w/Bearings - Rear	915-4457NA	915-4456NA	915-4457NA	915-4456NA
11	1	Parallel Arm Assy., Top, w/Bearings - Front	915-0422FA	915-5422FNA	915-0422FA	915-5422FNA
12	1	Parallel Arm Assy., Top, w/Bearings - Rear	915-5422RNA	915-0422RA	915-5422RNA	915-0422RA
13	1	Parallel Arm Assy., Bottom, w/Bearings - Front	915-6458NA	915-6455NA	915-6458NA	915-6455NA
14	1	Rubber Nose	24603-30	24603-30	24603-33	24603-33

Part Numbers of Items Identical on All Lift Models							
Item	Qty.	Description	Part No.	Item	Qty.	Description	Part No.
15	1	Block, Platform Stop	25778	100	2	Spring, Dampener, Retract	26963
16	1	Block, Platform Stop	900-0311	101	2	Bolt, 5/16-18 x 1 1/4" Low Socket Hd., Cap, Bk.	26327
17	4	Snap-Ring 5/16"	24570	102	2	Nut, 1/4-20, Hex	10057
18	4	O-Ring, 5/16" I.D. x 1/16" O.D.	26614	103	2	Screw, 1/2-20 x 1 1/2" Set Lock, Black	18663
19	2	Screw, Socket Head, 3/8-16 x 1/2" with 1/2" Shoulder	25929	104	2	Block, Guide-Platform Stow	28593
20	2	Grommet, I.D. .281"/FR.125" x GD.625"/T.375"	27725	105	2	Pin, Platform	975-2325
21	1	Pivot, Outer Barrier Lever	24550	106	10	Screw, 5/16-18 x 3/4", SHFS	25527
22	2	Rivet, Pop, 3/16" x .06"/.13"	11512	107	2	Spacer, UHMW, Vertical Channel	995-0320
23	1	Stud, Wing Head	28804	108	1	Position Sensor Collet, +/- 65°	31837
24	1	Retainer, Push-On	28805	109	2	Bracket, Mounting, Quiet Ride	915-0392
25	1	Grommet, ID .312", FR .250", GD .813" T500"	27726	110	2	Gas Spring, 12.2" Extended/8.3" Compressed	30550A
26	1	Pump Assembly	31038-12V	111	2	Bolt, 5/16-18 x 3/4", Carriage, Black	15858BK
27	10	Pin, Pivot, Parallel Arm	917-0403	112	3	Washer, 5/16", External Tooth	16368
28	10	Screw Drive #8 x 3/16" Long Black	31086BK	113	4	Nut, 5/16-18 Hex	10058
29	1	Receptacle, Clip-On	28803	114	2	Bracket, Inner Side Panel Guide	915-0703
30	1	Cylinder Assembly, Outer Barrier	915-3318NA	115	3	Bearing, UHMW, Flat, Thin, Black	916-5406
31	4	Snap-Ring 1/2"	13273	116	2	Pivot, Inner Rollstop Lever, Upper, Inside	916-5433
32	1	Latch Assembly w/Magnet, Outer Barrier	915-4343NA	117	4	Bearing, Flange - 3/4" x 1/2"	24027
33	2	Bolt, 5/16-18 x 1/2" Nylock, Hex	10012	118	1	Spacer, UHMW, Latch	950-7760
34	2	Washer, 5/16" Flat	10063	119	1	Spring, 1/4" x 1" Extension w/Hooks	10774
35	2	Bolt, 3/8-16 x 1/2" Nylock, Hex	29729	120	1	Tubing, Polyurethane - 1/4" I.D. x 1/2" O.D. x .81"	28630R000.81
36	1	Assembly, IB Occupied	33779A	121	2	Pin, Clevis, 3/8" x 3", Black	24932BK
37	1	Assembly, Rotary Position Sensor, +/- 65°	31844A	122	1	Pin, Platform Cylinder	915-0414
38	2	Bolt, 5/16-18 x 3/4" Nylock, Hex	29608	123	2	Ring, Snap, 1/2", External, .035" Thick	20946
39	1	Rivet, Pop, SD68 BSLF-3/16" - .38"/.50"	15244	124	1	Weldment, Tube-IB Fold Arm-Front	918-0625RNW
40	5	Rivet, Pop, 1/8" - .25"/.31"	13906	125	1	Lever, Outer Barrier Latch (Assembly #915-0345A)	915-0345
41	4	Bearing, Flange, 1/2" x 1/4"	24442	126	4	Nut, 3/8-16, Hex Lock Jamb	20926
42	1	Weldment, Latch, Lift/Tite, Rear	27013RW	127	1	Weldment, Tube-IB Fold Arm-Rear	918-0625FNA
43	1	Weldment, Latch, Lift/Tite, Front	27013FW	128	1	Pivot, Outer Barrier Cylinder, Front	990-0352
44	2	Screw, #4-40 x 5/8", Round Head	14810	129	2	Bolt, 5/16-18 x 1 1/2" Hex Head, Cap	10014
45	2	Screw, 1/4-20 x 3/4", Flat Head, Countersink	28986	130	2	Bolt, 7/16-20 x 1", SHCS, w/Nylock Patch	28474
46	1	Screw, #8-32 x 1", Mach. Truss Head	19912BK	131	2	Screw, #10-32 x 3/8", Flat Head, Hex	24537
47	2	Nut, #4-40 w/Lockwasher	19537	132	3	Adapter, Cylinder/Outer Barrier Lever-Out	990-0341
48	2	Nut, #8-32, Hex	10777	133	2	Bolt, 5/16-18 x 3/4", Button Head, Soc Cap Hd	24440
49	1	Bag, Vinyl - 7" x 10"	31174	134	7	Rivet, Pop, 3/16"- .25"/.38", Black	14993
50	2	Cylinder, 14.625"/23.146" Retracted	C1514-3-9407	135	1	Adapter, Cylinder/Outer Barrier Lever	915-0341N
51	1	Screw, 10-32 x 5/8" Pan Head Phillips, Thd. Cut	26714	136	1	Harness, Sensor - Cherry w/Plug	32645NA60
52	4	Washer, 1/4" External Tooth	83588	137	4	Screw, #4-40 x 1/4", Phillips Head	11480
53	1	Clamp, 5/8" Closed Insulated	84070	138	11	Clip, Cable, 1/2", Plastic	10092
54	1	Vertical Channel Weldment, Front, 42"	918-0442FNW	139	2	Nut, 5/16-18, Top Lock, Black	28324BK
55	1	Cover, Plastic, Parallel Arm, O/S Front-w/Tape	915-0704NA	140	2	Bolt, 1/4-20 x 1/2", Button Head Socket	15733
56	2	Tape, 1" Wide S-Face Foam, x 14" (1 Shown)	82015R014	141	2	Bolt, 3/8-16 x 2", Hex Head, Cap	10027
57	6	Tape, 1" Wide S-Face Foam, x 10" (3 Shown)	82015R010	142	2	Roller Assembly	1000-2395A
58	1	Cycle Counter-LCD Without Reset	30547	143	2	Roller Retainer	915-5353
59	2	Handrail Weldment	985-4618NW	144	2	Nut, 5/16-18 Hex, Lock	11387BK
60	8	Rivet, Push In - 8MM	30063	145	2	Pivot, Inner Rollstop Lever, Upper-Outer	916-5434
61	4	Bumper, Rubber, 31/32" O.D.	10950	146	6	Clamp, 1/4" I.D. - Nylon Loop	84396
62	2	Screw, #8-18 x 1/2", Self-Tap, Phillips	13583	147	3	Washer, #10, Flat, Auto BK	11541
63	2	Bulb, Light 20-W, Halogen Reflector	31060	148	18	Rivet, Snap, .201" x .256"/.295"	26362
64	2	Nut, #10-32 w/Lockwasher	18349	149	4	Bearing, Flange, 1" x 1/2-16 FDU08	28031
65	2	Screw, 1/4-20 x 3/4", Button Head, Socket Cap	81068-000	150	1	Cover, Plastic, Parallel Arm, O/S Rear-w/Tape	915-0702NA
66	6	Washer, 3/8" Flat, Black	10064	151	1	Weldment, Inner Rollstop Bracket, Front	985-0431FNW
67	4	Nut, 3/8-16, Hex Lock	13617	152	4	Spacer, Handrail, .375" Thk. UHMW	995-0301
68	2	Cotter Pin, 1/4" x 2", Black	15886BK	153	2	Washer, .516" I.D. x .75" O.D. x .250", Black	25336
69	2	Nut, 1/4-20 Hex Lock	10775	154	1	Cover, Plastic, Parallel Arm, I/S Rear-w/Tape	975-0703NA
70	1	Cover, Plastic, Parallel Arm, I/S Front-w/Tape	975-0701NA	155	1	Weldment, Inner Rollstop Bracket, Rear	985-0431RNW
71	2	Threshold Strip Switch	31221A	156	2	Spring, Torsion - .366 I.D. x .05"	30386
72	2	Bolt, Shoulder, 5/16" x 5/8"-1/4-20	31816	157	4	Spacer, IB Lock	30394
73	10	Washer, .516" I.D. x 1" O.D. x .0269", Auto BK	15328	158	2	Weldment, IB Lock	985-0630NW
74	10	Washer, 1/2" Flat Black Zinc	25346	159	4	Bearing, 3/8" I.D. x 1/2" Long	30405
75	2	Bolt, 3/8-16 x 1", Hex	10025	160	4	Roller, .65" O.D. x 1/2" Long	985-0631N
76	4	Weld Nut, #4-40 x 1/8" Barrel	32949	161	4	E-Clip, 3/8" Shaft	84383
77	2	Lever, Inner Rollstop	985-0404N	162	2	Spacer, Handrail, .230" Thk. UHMW	911-5301
78	4	Clamp, Insulate, 1 3/8" I.D.	29765	163	1	Adapter, Cylinder/Outer Barrier Lever	915-0344
79	2	Washer, 5/16" Lock, Black	10068	164	1	Nut, 1/2-20, Hex, Jam	83022
80	2	Gas Spring, Fold Arm	29186	165	3	Rivet, Pop, 3/16" - .13/.25" - SD64BS	11513
81	2	Ball Stud, Gas Spring, 13MM x 3/8-16 F. Thread	29185	166	2	Bearing, Flange, 3/8" x 1/4" - Plastic	24028
82	2	Washer, 3/8", Lock	10069	167	2	Light Assy, Vertical Channel	31062A
83	1	Spacer, Flap Cylinder .184" Thick	915-0416	168	1	Tab, Outer Barrier Assembly	915-6367NA
84	2	Clip, Safety, 3/8-16 x 3/4", Flange	12608	169	1	Washer, .4" I.D. x .810" O.D. x .12"	24536
85	6	Bolt, 3/8-16 X 3/4", Flange Button Head, Socket Cap	25171	170	4	Spacer, UHMW-.75" O.D. x .39" I.D. x .25"	30227
86	2	Bolt, 5/16-18 x 2", Flange, Hex	25130	171	1	Bolt, 3/8-16 x 3/4" Soc. Hd. Cap Screw Bk.	26259
87	2	Nut, 5/16-18, Serrated Flange	25131	172	1	Spacer, Flap Cylinder .700" Thick	915-0415
88	2	Arm, Slide Support (Replace In Pairs)	985-6612N	173	1	Guard, Outer Barrier Sensor Wires	915-0368N50Y
89	2	Slide, Platform, Rotate, Weldment	915-0606N	174	6	Washer, Finger, Disc Spring	27276
90	2	Spacer, .5" O.D. x .334" I.D. x 1.515" LG	25132	175	4	Screw, #10-32 x 1/2" Button Head, Socket, Cap	30375WP
91	10	Bearing, Flange, 3/4" x 1/4"	24012	176	1	Position Sensor Spacer	31092
92	2	Weldment, Platform Pivot Arm - 42"	918-2640NW-BKN	177	1	Slide, UHMW, Platform Rotate	31677
93	2	Pin, Pivot, Lower Arm	900-0413N	178	2	Bolt, 1/4-20 x 1", Hex	10001
94	6	Ring, Snap, 3/4", External	18657	179	2	Bumper, 1" x 1/2", UHMW	33654
95	2	Bearing, UHMW, Flat, 1.226"	205-1760	180	2	Washer, 3/8" Lock, Black	10069BK
96	2	Ball Stud, 10MM	21301	181	2	Screw, #10-32 x 1/2" FHDHXS, Black w/ Patch	17192P
97	1	Vertical Channel Weldment, Rear, 42"	918-0442RNW				
98	6	Plastic End Cap, 1" I.D. x 1 1/4" O.D.	20474				
99	2	Grip, Handle, Yellow	19074				

42" Floor-to-Ground Lift Exploded View

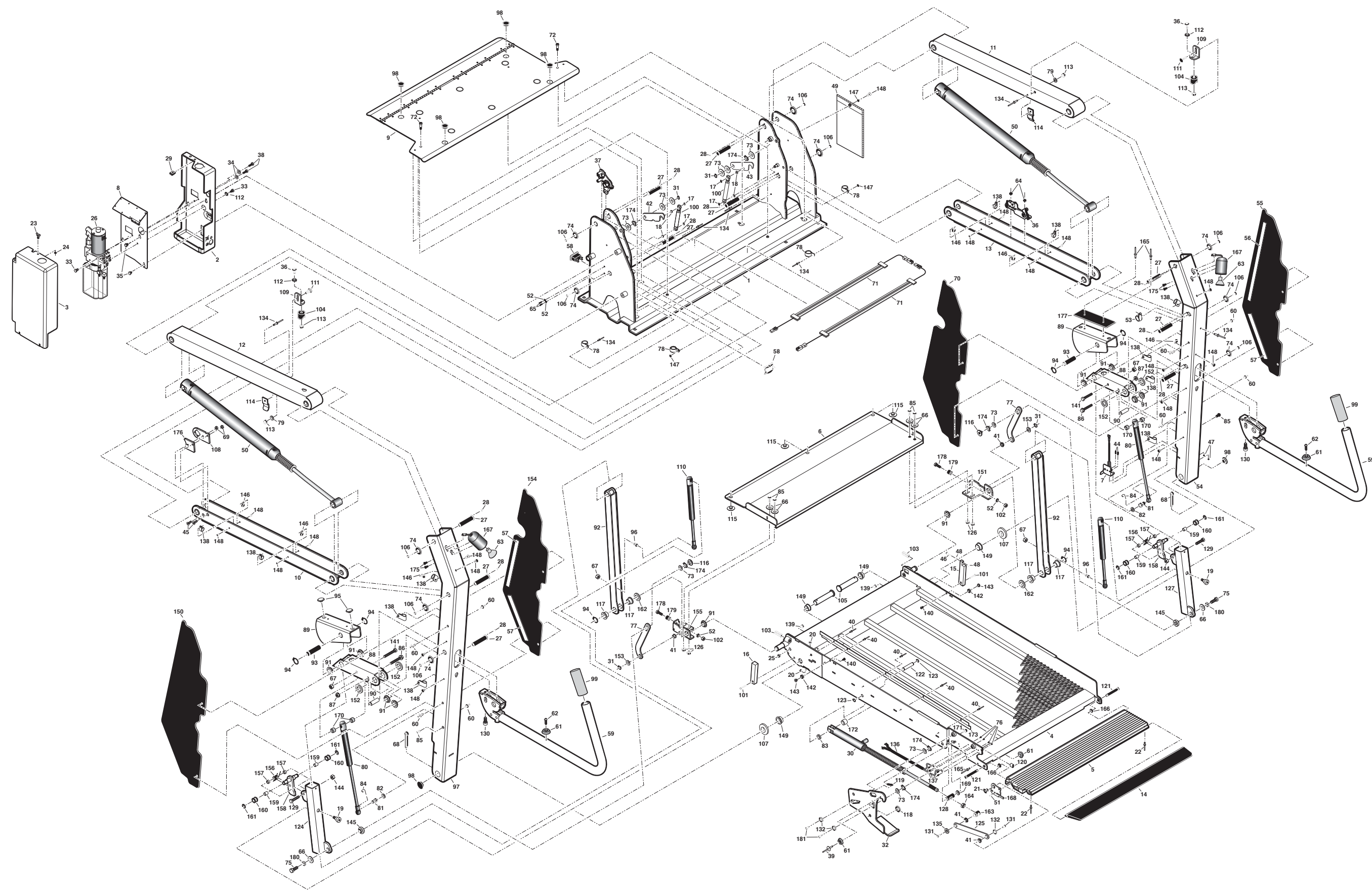

**Unfold for:
42" F-T-G Lift
Exploded
View**

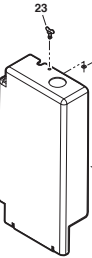

48" Floor-to-Ground Lift Repair Parts List

Part Numbers of Items Dedicated per Lift Model							
Item	Qty.	Description	NL916IB	NL916FIB	NL917IB	NL917FIB	NL919IB
1	1	Base Weldment	985R5148NW30	985F5148NW30	985R5148NW33	985F5148NW33	985R5148NW34
2	1	Cover, Pump, 2-Piece, Back-Bottom	915-0513RNPT	915-0513FNPT	915-0513RNPT	915-0513FNPT	915-0513FNPT
3	1	Cover, Pump, 2-Piece, Top-Front	915-0519RNPT	915-0519FNPT	915-0519RNPT	915-0519FNPT	915-0519FNPT
4	1	Platform Weldment	917-53050NW	917-53050NW	917-53350NW	917-53350NW	919-53450NW
5	1	Outer Barrier Assembly (Inc. 5, 10, 31, 165, 168)	915-0312NA30Y	915-0312NA30Y	915-0312NA33Y	915-0312NA33Y	915-0312NA34Y
6	1	Inner Roll Stop	915-0147NW30Y	915-0147NW30Y	915-0147NWY	915-0147NWY	915-0147NW34Y
7	1	Bridge Switch Assembly	31010FA	31010FA	31010FA	31010FA	31010FA
8	1	Plate, Backing/Mounting	985-2501RN	985-2501FN	985-2501RN	985-2501FN	985-2501RN
9	1	Base Cover Assy., (Traction Tape Not Shown)	975-5148CNA30Y	975-5148CNA30Y	975-5148CNA33Y	975-5148CNA33Y	975-5148CNA34Y
10	1	Rubber Nose	24603-30	24603-30	24603-33	24603-33	24603-34
11	1	Parallel Arm Assy., Top-Front w/Bearings	915-0437FA	985-5437FNA	915-0437FA	985-5437FNA	915-0437FA
12	1	Parallel Arm Assy., Top Rear w/Bearings	985-5437RNA	917-0436RNA	985-5437RNA	917-0436RNA	985-5437RNA
13	1	Parallel Arm Assy., Bottom Rear w/Bearings	917-4458RNA	917-4457NA	917-4458RNA	917-4457NA	917-4458RNA
14	1	Parallel Arm Assy., Bottom Front w/Bearings	917-6458FNA	975-6458NA	917-6458FNA	975-6458NA	917-6458FNA
Part Numbers of Items Identical on All Lift Models							
Item	Qty.	Description	Part No.	Item	Qty.	Description	Part No.
15	1	Block, Platform Stop - Bridging - Rear	25778	100	2	Spring, Dampener, Retract	26963
16	1	Block, Platform Stop - Bridging - Front	900-0311	101	2	Bolt, 5/16-18 x 1 1/4" Low Socket Hd., Cap, Bk.	26327
17	4	Snap-Ring 5/16"	24570	102	1	Nut, 1/2-20, Hex, Jam	83022
18	4	O-Ring, 5/16" I.D. x 1/16" O.D.	26614	103	2	Screw, 1/2-20 x 1 1/2" Set Lock, Black	18663
19	2	Screw, Socket Head, 3/8-16 x 1/2" with 1/2" Shoulder	25929	104	2	Block, Guide-Platform Stow	28593
20	2	Grommet, I.D. .281" / FR.125" x GD.625"/T.375"	27725	105	2	Pin, Platform	975-2325
21	1	Pivot, Outer Barrier Lever	24550	106	10	Screw, 5/16-18 x 3/4", SHFS	25527
22	1	Slide, UHMW, Platform Rotate	31677	107	2	Spacer, UHMW, Vertical Channel	995-0320
23	1	Stud, Wing Head	28804	108	1	Position Sensor Collet	31093
24	1	Retainer, Push-On	28805	109	2	Bracket, Mounting, Quiet Ride	915-0392
25	1	Grommet, ID .312", FR. .250", GD .813" T500"	27726	110	2	Gas Spring, 12.2" Extended/8.3" Compressed	30550A
26	1	Pump Assembly	31038-12V	111	2	Bolt, 5/16-18 x 3/4", Carriage, Black	15858BK
27	10	Pin, Pivot, Parallel Arm	917-0403	112	3	Washer, 5/16", External Tooth	16368
28	10	Screw Drive #8 x 3/16" Long Black	31086BK	113	4	Nut, 5/16-18 Hex	10058
29	1	Receptacle, Clip-On	28803	114	2	Bracket, Inner Side Panel Guide	915-0703
30	1	Cylinder Assembly, Outer Barrier	915-3318NA	115	3	Bearing, UHMW, Flat, Thin, Black	916-5406
31	1	Screw, 10-32 x 5/8" Pan Head Phillips, Thd. Cut	26714	116	2	Pivot, Inner Rollstop Lever, Upper, Inside	916-5433
32	1	Latch, Outer Barrier Assembly w/Magnet	915-4343NA	117	4	Bearing, Flange - 3/4" x 1/2"	24027
33	2	Bolt, 5/16-18 x 1/2" Nylock, Hex	10012	118	1	Spacer, UHMW, Latch	950-7760
34	2	Washer, 5/16" Flat	10063	119	1	Spring, 1/4" x 1" Extension w/Hooks	10774
35	2	Bolt, 3/8-16 x 1/2" Nylock, Hex	29729	120	1	Tubing, Polyurethane - 1/4" I.D. x 1/2" O.D. x .81"	28630R000.81
36	1	Assembly, IB Occupied	33779A	121	2	Pin, Clevis, 3/8" x 3", Black	24932BK
37	1	Assembly, Rotary Position Sensor	31094A	122	1	Pin, Platform Cylinder	915-0414
38	2	Bolt, 5/16-18 x 3/4" Nylock, Hex	29608	123	2	Ring, Snap, 1/2", External, .035" Thick	20946
39	1	Rivet, Pop, SD68 BSLF-3/16" - .38"/.50"	15244	124	1	Weldment, Tube-IB Fold Arm-Rear	985-0625RNW
40	5	Rivet, Pop, 1/8" - .25"/.31"	13906	125	1	Lever, Outer Barrier Latch (Assembly #915-0345A)	915-0345
41	4	Bearing, Flange, 1/2" x 1/4"	24442	126	4	Nut, 3/8-16, Hex	10059BK
42	1	Weldment, Latch, Lift/Tite, Rear	27013RW	127	1	Weldment, Tube-IB Fold Arm-Front	985-0625FNW
43	1	Weldment, Latch, Lift/Tite, Front	27013FW	128	1	Pivot, Outer Barrier Cylinder, Front	990-0352
44	2	Screw, #4-40 x 5/8", Round Head	14810	129	2	Bolt, 5/16-18 x 1 1/2" Hex Head, Cap	10014
45	2	Screw, 1/4-20 x 3/4", Flat Head, Countersink	28986	130	2	Bolt, 7/16-20 x 1", SHCS, w/Nylock Patch	28474
46	1	Screw, #8-32 x 1", Mach. Truss Head	19912BK	131	2	Screw, #10-32 x 3/8", Flat Head, Hex	24537
47	2	Nut, #4-40 w/Lockwasher	19537	132	3	Adapter, Cylinder/Outer Barrier Lever-Out	990-0341
48	2	Nut, #8-32, Hex	10777	133	2	Bolt, 5/16-18 x 3/4", Button Head, Soc Cap Hd	24440
49	1	Bag, Vinyl - 7" x 10"	31174	134	6	Rivet, Pop, 3/16" - .25"/.38", Black	14993
50	2	Cylinder, 14.625"/29.146" Retracted	C1514-3-9408	135	1	Adapter, Cylinder/Outer Barrier Lever	915-0341N
51	1	Adapter, Cylinder/Outer Barrier Lever	915-0344	136	1	Harness, Sensor - Cherry w/Plug	32645NA60
52	4	Washer, External Tooth	83588	137	4	Screw, #4-40 x 1/4" Phillips Head	11480
53	6	Washer, Finger, Disc Spring	27276	138	11	Clamp, Cable, 1/2", Plastic	10092
54	1	Vertical Channel Weldment, Front, 48"	955-0448FNW	139	2	Nut, 5/16-18, Top Lock, Black	28324BK
55	1	Cover, Plastic, Parallel Arm, O/S Front-w/Tape	915-0704NA	140	2	Bolt, 1/4-20 x 1/2", Button Head Socket	15733
56	2	Tape, 1" Wide S-Face Foam, x 14" (1 Shown)	82015R014	141	2	Bolt, 3/8-16 x 2", Hex Head, Cap	10027
57	6	Tape, 1" Wide S-Face Foam, x 10" (3 Shown)	82015R010	142	2	Roller Assembly	1000-2395A
58	1	Cycle Counter-LCD Without Reset	30547	143	2	Roller Retainer	915-5353
59	2	Handrail Weldment	985-4618NW	144	2	Nut, 5/16-18 Hex, Lock	11387BK
60	8	Rivet, Push In - 8MM	30063	145	2	Pivot, Inner Rollstop Lever, Upper-Outer	916-5434
61	4	Bumper, Rubber, 31/32" O.D.	10950	146	6	Clamp, 1/4" I.D. - Nylon Loop	84396
62	2	Screw, #8-18 x 1/2", Self-Tap, Phillips	13583	147	1	Clamp, 5/8" Closed Insulated	84070
63	2	Bulb, Light 20-W, Halogen Reflector	31060	148	18	Rivet, Snap, .201" x .256"/.295"	26362
64	2	Nut, #10-32 w/Lockwasher	18349	149	4	Bearing, Flange, 1" x 1/2-16FUDU08	28031
65	2	Screw, 1/4-20 x 3/4", Button Head, Socket Cap	81068-000	150	1	Cover, Plastic, Parallel Arm, O/S Rear-w/Tape	915-0702NA
66	2	Washer, 3/8" Flat, Black	10064	151	1	Weldment, Inner Rollstop Bracket, Front	985-0431FNW
67	4	Nut, 3/8-16, Hex Lock	13617	152	4	Spacer, Handrail, .375" Thk. UHMW	995-0301
68	2	Cotter Pin, 1/4" x 2", Black	15886BK	153	2	Washer, .516" I.D. x .75" O.D. x .250", Black	25336
69	2	Nut, 1/4-20 Hex Lock	10775	154	1	Cover, Plastic, Parallel Arm, I/S Rear-w/Tape	975-0703NA
70	1	Cover, Plastic, Parallel Arm, I/S Front-w/Tape	975-0701NA	155	1	Weldment, Inner Rollstop Bracket, Rear	985-0431RNW
71	2	Threshold Strip Switch	31221A	156	2	Spring, Torsion - .366 I.D. x .05"	30386
72	2	Bolt, Shoulder, 5/16 x 5/8-1/4"-20	31816	157	4	Spacer, IB Lock	30394
73	3	Washer, #10, Flat, Auto BK	11541	158	2	Weldment, IB Lock	985-0630NW
74	10	Washer, 1/2" Flat Black Zinc	25346	159	4	Bearing, 3/8" I.D. x 1/2" Long	30405
75	2	Bolt, 3/8-16 x 1", Hex	10025	160	4	Roller, .65" O.D. x 1/2" Long	985-0631N
76	6	Washer, 3/8", Lock, Black	10069BK	161	4	E-Clip, 3/8" Shaft	84383
77	2	Lever, Inner Rollstop	985-0404N	162	2	Spacer, Handrail, .230" Thk. UHMW	911-5301
78	4	Clamp, Insulate, 1 3/8" I.D.	29765	163	10	Washer, .516" I.D. x 1" O.D. x .269", Auto BK	15328
79	2	Washer, 5/16" Lock, Black	10068	164	3	Rivet, Pop, 3/16" - .13/25" - SD64BS	11513
80	2	Gas Spring, Fold Arm	29186	165	2	Rivet, Pop, 3/16" x .06"/.13	11512
81	2	Ball Stud, Gas Spring, 13MM x 3/8-16 F. Thread	29185	166	2	Bearing, Flange, 3/8" x 1/4" - Plastic	24028
82	2	Washer, 3/8", Lock	10069	167	2	Light Assy, Vertical Channel	31062A
83	1	Spacer, Flap Cylinder .184" Thick	915-0416	168	1	Tab, Outer Barrier - Assembly	915-6367NA
84	2	Clip, Safety, 3/8-16 x 3/4" Flange	12608	169	1	Washer, 4" I.D. x 810" O.D. x .12"	24536
85	6	Bolt, 3/8-16 x 3/4", Flange Button Head, Socket Cap	25171	170	4	Spacer, UHMW-.75" O.D. x .39" I.D. x .25"	30227
86	2	Bolt, 5/16-18 x 2", Flange, Hex	25130	171	1	Bolt, 3/8-16 x 3/4" Soc. Hd. Cap Screw Bk.	28259
87	2	Nut, 5/16-18, Serrated Flange	25131	172	1	Spacer, Flap Cylinder .700" Thick	915-0415
88	2	Arm, Slide Support (Replace In Pairs)	985-6612N	173	1	Guard, Outer Barrier Sensor Wires	915-0368N50Y
89	2	Slide, Platform, Rotate	985-0606N	174	4	Snap-Ring 1/2"	13273
90	2	Spacer, .5" O.D. x .334" I.D. x 1.515" LG	25132	175	4	Screw, #10-32 x 1/2" Button Head, Socket, Cap	30375WP
91	10	Bearing, Flange, 3/4" x 1/4"	24012	176	1	Position Sensor Spacer	31092
92	2	Weldment, Platform Pivot Arm - 48"	985-0640NW-BKN	177	4	Weld Nut, #4-40 x 1/8" Barrel	32949
93	2	Pin, Pivot, Lower Arm	900-0413N	178	2	Bolt, 1/4-20 x 1", Hex	10001
94	6	Ring, Snap, 3/4", External	18657	179	2	Bumper, 1" x 5/8", UHMW	33603
95	2	Bearing, UHMW, Flat, 1.226"	205-1760	180	2	Screw, #10-32 x 1/2" FHDHXS, Black w/Patch	17192P
96	2	Ball Stud, 10MM	21301	181	2	Nut, 1/4-20, Hex	10057
97	1	Vertical Channel Weldment, Rear, 48"	955-0448RNW				
98	6	Plastic End Cap, 1" I.D. x 1 1/4" O.D.	20474				
99	2	Grip, Handle, Yellow	19074				

Braun® Limited Warranty

WARRANTY COVERAGE AND WARRANTY COVERAGE TIME PERIODS

The Braun Corporation ("Braun") warranty covers certain parts of this wheelchair lift for three (3) years or 10,000 cycles and the cost of labor to repair or replace those parts for one (1) year or 3,000 cycles. If The Braun Corporation receives the warranty registration card within 20 days after the lift is put into service, the warranty labor coverage will increase from one (1) year or 3,000 cycles to three (3) years or 10,000 cycles. In addition, providing the warranty registration card is returned as noted above, the following lift's power train parts are warrantied for five (5) years or 15,000 cycles: Cable, Cylinder, Flow Control, Gear Box, Motor, Pump, Hydraulic Hose and Fittings. This limited warranty covers substantial defects in materials and workmanship of the lift, provided that the lift is operated and maintained properly and in conformity with the owner's manual. The warranty period begins on the date that the product is delivered to the first retail purchaser by an independent, authorized dealer of Braun, or, if the dealer places the product into any type of service prior to retail sale, on the date the dealer first places the product in such service. This limited warranty applies only to the first purchaser. It may not be transferred.

WHAT BRAUN WILL DO TO CORRECT PROBLEMS

In the event that a substantial defect in material or workmanship, attributable to Braun, is found to exist during the first year of warranty coverage, it will be repaired or replaced, at Braun's option, without charge for parts or labor to the owner, in accordance with the terms, conditions and limitations of this limited warranty. If the substantial defect in material or workmanship, attributable to Braun, is found to exist during the second or third year of warranty coverage, it will be repaired or replaced, at Braun's option, without charge to the owner for parts, only, in accordance with the terms, conditions and limitations of this limited warranty. Providing the warranty card is returned within 20 days as outlined above, the labor warranty period will be extended by two years of coverage in accordance with the terms, conditions, and limitations of this limited warranty. In addition, if a substantial defect in material or workmanship, attributable to Braun, is found to exist during the fourth or fifth year of warranty coverage to the following lift's power train parts: Cable, Cylinder, Flow Control, Gear Box, Motor, Pump, Hydraulic Hose and Fittings, it will be repaired or replaced, at Braun's option, without charge to the owner for parts, only, in accordance with the terms, conditions and limitations of this limited warranty. The cost of labor for repair or replacement at any time after the warranty coverage detailed above is the sole responsibility of the owner.

Braun's obligation to repair or replace defective materials or workmanship is the sole obligation of Braun under this limited warranty. Braun reserves the right to use new or remanufactured parts of similar quality to complete any work, and to make parts and design changes from time to time without notice to anyone. Braun reserves the right to make changes in the design or material of its products without incurring any obligation to incorporate such changes in any previously manufactured product. Braun makes no warranty as to the future performance of this product, and this limited warranty is not intended to extend to the future performance of the product. In addition, the owner's obligation to notify Braun, or one of its authorized, independent dealers, of a claimed defect does not modify any obligation placed on the owner to contact Braun directly when attempting to pursue remedies under state or federal law.

LIMITATIONS, EXCLUSIONS AND DISCLAIMER OF IMPLIED WARRANTIES

ANY IMPLIED WARRANTY THAT IS FOUND TO ARISE BY WAY OF STATE OR FEDERAL LAW, INCLUDING ANY IMPLIED WARRANTY OF MERCHANTABILITY OR ANY IMPLIED WARRANTY OF FITNESS, IS LIMITED IN DURATION TO THE TERMS OF THIS LIMITED WARRANTY AND IS LIMITED IN SCOPE OF COVERAGE TO THE SCOPE OF COVERAGE OF THIS LIMITED WARRANTY. Braun disclaims any express or implied warranty, including any implied warranty of fitness or merchantability, on items excluded from coverage as set forth in this limited warranty. Braun makes no warranty of any nature beyond that contained in this limited warranty. No one has authority to enlarge, amend or modify this limited warranty, and Braun does not authorize anyone to create any other obligation for it regarding this product. Braun is not responsible for any representation, promise or warranty made by any independent dealer or other person beyond what is expressly stated in this limited warranty. Any selling or servicing dealer is not Braun's agent, but an independent entity.

Braun® Limited Warranty

BRAUN SHALL NOT BE LIABLE FOR ANY INCIDENTAL OR CONSEQUENTIAL DAMAGES THAT MAY RESULT FROM BREACH OF THIS LIMITED WARRANTY OR ANY IMPLIED WARRANTY. THIS EXCLUSION OF CONSEQUENTIAL AND INCIDENTAL DAMAGES SHALL BE INDEPENDENT OF ANY FAILURE OF THE ESSENTIAL PURPOSE OF ANY WARRANTY, AND THIS EXCLUSION SHALL SURVIVE ANY DETERMINATION THAT THIS LIMITED WARRANTY OR ANY IMPLIED WARRANTY HAS FAILED OF ITS ESSENTIAL PURPOSE. This warranty does not cover, and in no event shall Braun be liable for towing charges, travel, lodging, or any other expense incurred due to the loss of use of the product or other reason.

Some states do not allow limitations on how long an implied warranty lasts, or the exclusion or limitation of incidental or consequential damages, so the above limitations or exclusions may not apply to you.

HOW TO GET SERVICE

To obtain warranty service the owner must do all of the following:

1. Notify an authorized service center, of the claimed defect attributable to Braun, within the warranty coverage period designated above
2. Provide the notification mentioned in (1), above, within ten (10) days of when the owner discovered, or should have discovered, the claimed defect
3. Promptly schedule an appointment with and take the product to an authorized service center for service.
4. Pay any transportation costs and all expenses associated with obtaining warranty service.

Since Braun does not control the scheduling of service work at the independent dealerships you may encounter some delay in scheduling or completion of work. If you need assistance you may contact Braun, at 631 West 11th Street, Winamac, Indiana 46996; 1-800-THE-LIFT, (843-5438).

If two (2) or more service attempts have been made to correct any covered defect that you believe impairs the value, use or safety of the product, or if it has taken longer than thirty (30) days for repairs to be completed, you must, to the extent permitted by law, notify Braun directly, in writing, at the above address, of the unsuccessful repair(s) of the alleged defect(s) so that Braun can become directly involved in providing service pursuant to the terms of this limited warranty.

WHAT IS NOT COVERED

This Limited Warranty does not cover any of the following: defects in materials, components or parts of the product not attributable to Braun, any material, component or part of the product that is warranted by another entity (Note: the written warranty provided by the manufacturer of the material, component or part is the direct responsibility of that manufacturer); items that are added or changed after the product leaves Braun's possession; additional items installed at any dealership, or other place of business, or by any other party, other than Braun; normal wear, tear, usage, maintenance, service, periodic adjustments, the effects of condensation or moisture from condensation; mold or any damage caused by mold; imperfections that do not affect the product for its intended purpose; items that are working as designed but that you are unhappy with; problems related to mis-operation, misuse, mishandling, neglect or abuse, including failure to maintain the product in accordance with the owner's manual, or other routine maintenance such as inspections, lubricating, adjustments, tightening of screws, sealing, wheel alignments or rotating tires; damage due to accident or collision, including any acts of weather or damage or corrosion due to the environment; theft, vandalism, fire, or other intervening acts not attributable to Braun; damage resulting from tire wear or tire failure; defacing, scratches, dents or chips on any interior or exterior surface of the product, including those caused by rocks or other road hazards, damage caused by off road use, overloading or alteration of the product, or any of its components or parts.

Defects and/or damage to interior and exterior surfaces and other appearance items may occur at the factory or when the product is in transit. These items are usually detected and corrected at the factory or by

Braun® Limited Warranty

a dealer prior to delivery to the purchaser. You must inspect the product for this type of damage when you take delivery. If you find any such defect or damage you must notify the selling dealer, or Braun, at the time of delivery to have these items covered by this limited warranty and to have work performed on the items at no cost to you as provided by this limited warranty.

EVENTS DISCHARGING BRAUN FROM OBLIGATION UNDER WARRANTY

The following shall completely discharge Braun from any express or implied warranty obligation to repair or replace anything and void this warranty: misuse, neglect, collision, accidents, failure to provide routine maintenance (See Owner's Manual), unauthorized alteration, off road use, Acts of Nature, damage from weather or the environment, theft, vandalism, tampering, fire, explosions, overloading the product and odometer tampering.

LEGAL REMEDIES

Any action to enforce any portion of this limited warranty, or any implied warranty, must be commenced within six (6) months after expiration of the warranty coverage period designated above or the action will be barred because of the passage of time. Any performance of repairs shall not suspend this limitation period from expiring. Any performance of repairs after the warranty coverage period has expired, or performance of repairs regarding any thing excluded from coverage under this limited warranty shall be considered "good will" repairs, and they will not alter the terms of this limited warranty, or extend the warranty coverage period or the filing limitation period in this paragraph. In addition, since it is reasonable to expect that the product will need some service during the warranty period; this warranty does not extend to future performance. It only sets forth what Braun will do and does not guarantee anything about the product for any time period. Nothing in this warranty, or any action of Braun, or any agent of Braun, shall be interpreted as an extension of any warranty period or the filing limitation period in this paragraph. Some states do not allow a reduction in the statute of limitations, so this reduction may not apply to you.

WARRANTY REGISTRATION and MISCELLANEOUS

Your warranty registration records should be completed and delivered to the appropriate companies, including the Braun Delivery Checklist & Warranty form. That form must be returned to Braun within twenty (20) days of purchase. The Braun warranty will not be registered unless this warranty registration is completed and received by Braun. Failure to file this warranty registration with Braun will not affect your rights under this limited warranty as long as you can present proof of purchase, but it can cause delays in obtaining the benefits of this limited warranty, and it changes the start date of the warranty to the date of final assembly of the product by Braun.

Braun agrees to repair or replace any of its factory installed parts found to have substantial defects within the appropriate warranty period designated above, provided that the repair is authorized by Braun and carried out by an authorized service center (a Braun labor schedule determines the cost allowance for repairs). Braun will not honor any warranty claim for repairs or replacement of parts unless the claim is submitted with the appropriate paperwork, and the work is completed by an independent, factory authorized service center. The appropriate paperwork can be obtained by written or phone contact with Braun at the contact information in this warranty.

Braun reserves the right to designate where any warranty work can be performed. Braun also reserves the right to examine any defective workmanship or part prior to giving any authorization for warranty work. Braun's return authorization procedure must be adhered to in order to process any warranty claims.

THIS WARRANTY GIVES YOU SPECIFIC LEGAL RIGHTS. YOU MAY ALSO HAVE OTHER RIGHTS THAT VARY FROM STATE TO STATE.

"Providing Access to the World"®



**Over 300 Braun
Dealers Worldwide**

 **THE BRAUN
CORPORATION®**
"Providing Access to the World"



International Corporate Hdqrs: P.O. Box 310 Winamac, IN 46996 USA
1-800-THE LIFT® (574) 946-6153 FAX: (574) 946-4670

Service Manual for:

NL

Millennium

Series

Public Use Wheelchair Lifts

Series C5

Patent #5,261,779

Patent #6,599,079

Patent #7,422,408

Patent #6,065,924

Patent #6,692,217

Patent #7,444,995

Patent #6,238,169

Patent #6,739,824

Patent #7,509,187

Patent #6,464,447

Patent #7,306,422

Patent #7,530,226

Patents Pending



"Providing Access to the World"®

International Corporate Hdqrs: P.O. Box 310 Winamac, IN 46996 USA
1-800-THE LIFT® (574) 946-6153 FAX: (574) 946-4670



33995 Rev. C
December 2009

BRAUN
MILLENNIUM
SERIES

Braun NL Series

All illustrations, descriptions and specifications in this manual are based on the latest product information available at the time of publication. The Braun Corporation reserves the right to make changes at any time without notice.

© The Braun Corporation